

i-STAT 1

Important!

i-STAT 1 System Manual Update

As of November 2022, the current i-STAT 1 System Manual has been updated. Please **ADD** and **DELETE** the sheets as listed below. Once the updates have been completed these instructions may be discarded. If you have any questions about these instructions, please contact your i-STAT Support Services provider.

A ► ADD SHEET

◀ D DESTROY SHEET

| <u>Item</u> | <u>Art#</u> |
|--|------------------------|
| A ► i-STAT 1 System Manual Configuration Sheet | 714419-00BK |
| ◀ D i-STAT 1 System Manual Configuration Sheet | 714419-00BJ (or lower) |
| <u>System Components Tab</u> | |
| A ► i-STAT 1 System Manual Section 1: Introduction | 714363-00AE |
| ◀ D i-STAT 1 System Manual Section 1: Introduction | 714363-00AD (or lower) |
| A ► i-STAT 1 System Manual Section 2: i-STAT 1 Analyzer | 714364-00X |
| ◀ D i-STAT 1 System Manual Section 2: i-STAT 1 Analyzer | 714364-00W (or lower) |
| <u>Care and Software Updates Tab</u> | |
| A ► Technical Bulletin: Network Options for Updating the i-STAT 1 Analyzer Using www.globalpointofcare.abbott | 731336-00J |
| ◀ D Technical Bulletin: Network Options for Updating the i-STAT 1 Analyzer Using www.globalpointofcare.abbott | 731336-00I (or lower) |
| <u>Technical Bulletins Tab</u> | |
| ◀ D New Ultralife 9-Volt Lithium Battery For Use With the i-STAT System | 730271-00B (or lower) |
| A ► Technical Bulletin: Support Services | 716144-00AM |
| ◀ D Technical Bulletin: Support Services | 716144-00AL (or lower) |

END

© 2022 Abbott Point of Care Inc., 100 and 200 Abbott Park Road, Abbott Park, IL 60064 • USA
i-STAT is a trademark of Abbott.



i-STAT 1 SYSTEM MANUAL CONFIGURATION

Please ensure that the contents of your System Manual are complete and up to date. In the event that your System Manual does not contain the current configuration, it is recommended that you contact your i-STAT support provider.

As of November 2022, your i-STAT 1 System Manual should be configured with the contents as listed below and in the order shown.

| ITEM | Art # |
|--|--------------|
| Cover Sheet | 714336-00R |
| Configuration Sheet | 714419-00BK |
| Table of Contents..... | 714362-00AH |
| Section 1..... | 714363-00AE |
| Section 2..... | 714364-00X |
| Section 3..... | 714365-00J |
| Section 4..... | 714366-00D |
| Section 5..... | 714367-00I |
| Section 6..... | 714368-00N |
| Section 7..... | 714369-00O |
| Section 8..... | 714370-00G |
| Section 9..... | 714371-00J |
| Section 10..... | 714372-00R |
| Section 11..... | 714373-00K |
| Section 12..... | 714374-00P |
| Technical Bulletin: Sample Type Customization on the i-STAT 1 Analyzer..... | 765893-00B |
| Section 13..... | 714375-00G |
| Section 14..... | 714376-00Z |
| Section 15..... | 714377-00U |
| Section 16..... | 714378-00J |
| Section 17..... | 714379-00K |
| Section 18..... | 714380-00R |
| Technical Bulletin: Instructions for Updating i-STAT 1 Analyzer Software Using www.globalpointofcare.abbott | 731335-00I |
| Technical Bulletin: Network Options for Updating the i-STAT 1 Analyzer Using www.globalpointofcare.abbott | 731336-00J |
| Section 19..... | 714381-00L |
| Technical Bulletin: Analyzer Coded Messages..... | 714260-00Z |
| Section 20..... | 714382-00F |

CTI Sheets

Introduction714258-00V

Technical Bulletins

Instructions for Restoring Analyzers that Produce *** for
Hematocrit and Quality Check Code 23.....721215-00D
K₂EDTA and K₃EDTA Customization for
Hematocrit on the i-STAT System.....716240-00F
ACT Test Result Options: Prewarmed vs. Non-Prewarmed
Result Calibration Modes for the i-STAT 1 Analyzer.....715617-00F
i-STAT Celite ACT and i-STAT Kaolin ACT
Heparin Linearity Procedure.....714547-00H
Support Services..... 716144-00AM
The i-STAT System and Waived Status.....731678-00E
The Presence of Latex Rubber in i-STAT
System Components.....721296-00G

INTRODUCTION 1

This Manual

This manual describes the i-STAT 1 Analyzer and the Data Manager software. Related sections are grouped behind tabs.

Note: Not all products are available in all regions.

Intended Use

The i-STAT 1 Analyzer is intended for use in the *in vitro* quantification of various analytes in whole blood or plasma in point of care or clinical laboratory settings. Analyzers and cartridges should be used by healthcare professionals trained and certified to use the system and should be used according to the facility's policies and procedures.

The i-STAT System is for *in vitro* diagnostics use. Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.

FDA Test Categorization

With the i-STAT 1 System, the FDA has categorized the tests included on the i-STAT G and Crea cartridges as waived when testing is performed using venous whole blood samples collected in lithium heparin evacuated tubes. Other venous whole blood samples, capillary and/or arterial samples tested using these same cartridges on the i-STAT 1 System are categorized by the FDA as moderate complexity.

For waived testing, laboratories are required to follow the manufacturer's requirements for the testing. They may elect to perform additional quality control testing (such as the QC required for a moderate complexity test) but this does not change the FDA categorization of the test as waived or release the laboratory's responsibility to follow the manufacturer's instructions for it as a waived test.


Other testing performed with the i-STAT 1 System (other than the testing performed using the aforementioned cartridges) is FDA categorized as "moderate complexity".

Overview of the i-STAT System

The i-STAT System incorporates a comprehensive group of components needed to perform blood analysis at the point of care. A portable handheld analyzer, a cartridge with the required tests, and 2-3 drops of blood will allow the caregiver to view quantitative test results for blood gas, chemistry and coagulation tests in approximately two minutes.

Portable printers and infrared communication devices allow all patient information obtained at the bedside to be printed on demand and transmitted to centralized information systems for record keeping and billing.

The Data Manager provides system management tools including real-time monitoring of testing and operator competency.

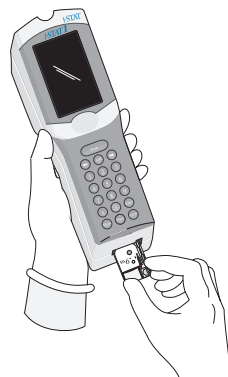

Components

The i-STAT System consists of:

- ✧ i-STAT Cartridges
- ✧ i-STAT 1 Analyzer
- ✧ Portable Printer
- ✧ Quality Assurance Materials
 - Electronic Simulator
 - Control Solutions
 - Calibration Verification Set (for cartridges)
- ✧ Data Management System
 - i-STAT 1 Downloader
 - i-STAT 1 Downloader/Recharger
 - Data Manager
- ✧ LIS/HIS Interface Software

Selection of Components

The selection of system components is dependent on factors unique to each facility such as:

- ✧ Types of tests to be performed
- ✧ Number of testing sites
- ✧ Number of tests per site
- ✧ System administration requirements

Summary of the Procedure








The cartridge test cycle is initiated by selecting i-STAT Cartridge from the Test menu or Quality Tests from the Administration menu. To perform cartridge testing, the operator fills a cartridge with sample, seals the cartridge with its snap or slide closure, and inserts the cartridge into the analyzer. The unit-use cartridge contains all applicable components to perform one or more tests including: calibrating solution, sample handling system, sensors and reagents. The analyzer automatically controls all steps in the testing cycle, which may include: fluid movement, reagent mixing, calibration and thermal control. Quality checks are performed continuously throughout the test cycle. Operator and patient IDs and patient chart information can be entered. When the test cycle is completed, results are displayed and the test record is stored.



















Data Management Test records can be transmitted to the Data Manager where they can be printed and/or transmitted to the Laboratory Information System or Hospital Information System. An optional portable printer enables the operator to print results at the point of care.









Interfacing The Data Manager can be interfaced to a Laboratory Information System (LIS) or Hospital Information System (HIS) to automate billing and patient record keeping.

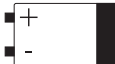


Note Regarding System Reliability The i-STAT System automatically runs a comprehensive set of quality checks of analyzer and cartridge performance each time a sample is tested. This internal quality system will suppress results if the analyzer or cartridge does not meet certain internal specifications (see Theory section in System Manual for detailed information). To minimize the probability of delivering a result with medically significant error the internal specifications are very stringent. It is typical for the system to suppress a very small percentage of results in normal operation given the stringency of these specifications. If however the analyzer or cartridges have been compromised, results may be persistently suppressed, and one or the other must be replaced to restore normal operating conditions. **Where unavailability of results while awaiting replacement of analyzers or cartridges is unacceptable, Abbott Point of Care Inc. recommends maintaining both a backup i-STAT System analyzer and cartridges from an alternate lot number.**

Symbols Symbols can be helpful in reducing the necessity for translating important information into multiple languages, particularly where space is limited. The following symbols may be found on components of the i-STAT System.
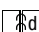
| Symbol | Definition |
|---|--|
|  | Attention: See instructions for use. |
|  | Caution: Risk of electrical shock. |
|  | Laser radiation hazard symbol. |
|  | Biological Risks. |
|  | Temperature limitations. The upper and lower limits for storage are adjacent to upper and lower arms. |
|  | Upper limit of temperature. The upper limit for storage is adjacent to the upper arm. |
|  | Use by or expiration date. An expiration date expressed as YYYY-MM-DD means the last day the product can be used. An expiration date expressed as YYYY-MM means the product cannot be used past the last day of the month specified. |

| Symbol | Definition |
|---|--|
|  | Manufacturer's lot number or batch code. The lot number or batch will appear adjacent to this symbol. |
|  | Catalog number, list number, or reference number. The number adjacent to this symbol is used to reorder the product. |
|  | Serial number. The serial number will appear adjacent to this symbol. |
| MN or # | Model number. The model number will appear adjacent to this symbol. |
|  | Date of manufacture. |
|  | Manufacturer. |
|  | <i>In vitro</i> diagnostic medical device. |
|  | For prescription use only. |
|  | Authorized Representative for Regulatory Affairs in the European Community. |
|  | Importer in the European Community. |
|  | A mark that indicates conformity to the legal requirements of the appropriate European Union (EU) Directive(s) with respect to safety, health, environment and consumer protection. |
|  | Contains sufficient for < n > tests. |
|  | Direct Current (DC). |
|  | Alternating Current (AC). |
|  | Class II Construction. |
|  | Consult instructions for use or see System Manual for instructions. |
|  | Control. |
|  | Signifies that the product bearing the ETL Listed mark complies with both U.S. and Canadian product safety standards: UL 61010-1: 3rd Ed.; Am. 1 CAN/CSA C22.2 No. 61010-1-12 3rd Ed. (R2017) +U1;U2 |
|  | i/immuno: Cartridges bearing this symbol must be run on i-STAT analyzers that also bear this symbol. |

| Symbol | Definition |
|---|--|
|  | Battery: i-STAT 1 Analyzer low battery icon (flashes on lower left side of display screen). |
|  | <p>Note concerning batteries: The following information is applicable to EEA (European Economic Area) countries: The directive 2006/66/EC requires separate collection of spent batteries. You are requested to dispose those batteries referred to on page 2-3 in accordance with local regulations.</p> <p>This product also contains a separate internal lithium battery that is not intended to be replaced by the user. See page 2-4 under "Additional Power" for more information.</p> |
|  | Separate waste collection for this electrical/electronic item indicated; Equipment manufactured / put on the market after 13 August 2005; Indicates compliance with Article 10(3) of Directive 2002/96/EC (WEEE) for the European Union (EU). |
| BODYYYY-MM-DD | Born On Date: the label BODYYYY-MM-DD defines year, month and day of manufacture. |
|  | Do not reuse. |
|  | This symbol is used for compliance with the China RoHS regulation (s). It indicates in years the Environmentally Friendly Use Period (EFUP) for the labeled electronic medical device product. |
| << >> | As the Martel Printer is incapable of printing the ↑ or ↓ symbols, this symbol appears on the Martel printout next to results which are outside the action range limits. |
| 14  | 14 days room temperature storage at 18-30 °C |
| 2  | 2 months room temperature storage at 18-30 °C |
| FCC | Signifies that the product bearing the Federal Communications Commission (FCC) logo complies with the specific requirements set forth by the FCC under Rules and Regulations, Title 47, Part 15 Subpart B, for Class A devices. |
|  | The near-patient testing symbol illustrates that a device can only be used in a near-patient setting by a health care worker, professional or trainee. |

| Symbol | The following symbols are used on the i-STAT 1 keypad. |
|---|--|
| SCAN | Key used to scan information into the analyzer. |
| ABC | Key used to enter letters. |
|  | Key used to enter information. |
| MENU | Key used to access the analyzer's menu. |
|  | Key used to print a test record. |
|  | Key used to turn the analyzer off and on. |

| Acronym | The following acronyms are listed in the i-STAT 1 System Manual |
|------------|---|
| CTI | Cartridge and Test Information |
| IFU | Instructions for Use |

| Symbol | The following symbols are used on i-STAT Value Assignment Sheets |
|-----------|--|
| \bar{x} | Mean |
| R | Range |

| Symbol | TEST |
|----------------------------|--|
| Na | Sodium |
| K | Potassium |
| Cl | Chloride |
| Glu | Glucose |
| Lac | Lactate |
| Crea | Creatinine |
| pH | pH |
| PCO₂ | Partial pressure of carbon dioxide |
| PO₂ | Partial pressure of oxygen |
| iCa | Ionized Calcium |
| BUN/UREA | Urea nitrogen/Urea |
| Hct | Hematocrit |
| ACTc Celite ACT | Activated Clotting Time with Celite® activator |

| Symbol | TEST |
|----------------------------|---|
| ACTk Kaolin ACT | Activated Clotting Time with Kaolin activator |
| PT/INR | Prothrombin Time / International Normalized Ratio |
| Hb | Hemoglobin |
| TCO₂ | Total carbon dioxide concentration |
| HCO₃ | Bicarbonate |
| BE (b&ecf) | Base excess (b for blood, ecf for extra cellular fluid) |
| AnGap | Anion Gap |
| sO₂ | Oxygen saturation |
| cTnI | Cardiac Troponin I |
| CK-MB | Creatine Kinase MB Isoenzyme |
| BNP | B-type Natriuretic Peptide |
| Total β-hCG | Total Beta-Human Chorionic Gonadotropin |

Warranty

Abbott Point of Care Inc. warrants this medical product (excluding disposable or consumable supplies) against defects in materials and workmanship for one year from the date of shipment. If Abbott Point of Care Inc. receives notice of such defects during the warranty period, Abbott Point of Care Inc. shall, at its option, either repair or replace products which prove to be defective. With respect to software or firmware, if Abbott Point of Care Inc. receives notice of defects in these products during the warranty period, Abbott Point of Care Inc. shall repair or replace software media and firmware which does not execute their programming instructions due to such defects. Abbott Point of Care Inc. does not warrant that the operating of the software, firmware or hardware shall be uninterrupted or error free. If Abbott Point of Care Inc. is unable, within a reasonable time, to repair or replace any product to a condition as warranted, Buyer shall be entitled to a refund of the purchase price upon return of the product to Abbott Point of Care Inc.

The warranty for the 9V rechargeable battery remains in effect for one year following the Born on Date (BODYYYY-MM-DD) noted on the battery's label (pictured below).



Note: Warranty rights may vary from state to state, province to province and country to country.

Limitations of Warranty

The foregoing warranty shall not apply to defects resulting from:

1. Improper or inadequate maintenance by Buyer or an unauthorized person,
2. Using accessories and/or consumables that are not approved by Abbott Point of Care Inc.,
3. Buyer-supplied software or interfacing,
4. Unauthorized repairs, modifications, misuse, or damage caused by disposable batteries, or rechargeable batteries not supplied by Abbott Point of Care Inc.
5. Operating outside of the environmental specifications of the product, or
6. Improper site preparation or maintenance.

THE WARRANTY SET FORTH ABOVE IS EXCLUSIVE AND NO OTHER WARRANTY, WHETHER WRITTEN OR ORAL, IS EXPRESSED OR IMPLIED. ABBOTT SPECIFICALLY DISCLAIMS THE IMPLIED WARRANTIES OR MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

© 2022 Abbott. All rights reserved. All trademarks referenced are trademarks of either the Abbott group of companies or their respective owners.

INTRODUCTION

The i-STAT 1 Analyzer is used in conjunction with i-STAT cartridges for the simultaneous quantitative determination of specific analytes in whole blood.

For information on the analytes that can be measured using i-STAT Cartridges, refer to the Cartridge and Test Information (CTI) sheets or Instructions for Use (IFU), located at www.globalpointofcare.abbott.

BEFORE YOU USE THE ANALYZER

- Install Batteries** See the Care of the Analyzer section in this manual for the procedure to install the disposable batteries. If a rechargeable battery is to be used, the disposable batteries can be used while the rechargeable battery pack is charged in the Downloader/Recharger. Charge rechargeable batteries fully before use. See the i-STAT 1 Downloader section for this procedure. When using a rechargeable battery, store the disposable battery carrier for possible future use.
- Check Date and Time** Press the On/Off key and check that the date and time at the top of the display are correct. To change the date and time, see Administration Menu in this section.
- Check Software** **Caution:** New analyzers or analyzers that have been repaired and returned or replaced will have standard CLEW and application software. If a different CLEW and/or application software is in use in your facility, it must be installed in new, repaired or replaced analyzers before they are put into use. Check the Analyzer Status page for the installed CLEW and application software. For steps to verify the software, see *Procedures for Customization using the Analyzer Keypad* in Section 9 of this manual. See under “Standardization and Calibration” in section 3 of this manual for an explanation of CLEW.
- Customization** Analyzers can be customized for many site-specific testing requirements. See Section 9 in this manual for a list of customizable parameters and their default values. To change the customization profile via the analyzer keypad see “Customization” under “Administration” in this section of the manual. To change the customization profile via i-STAT/DE, see the “i-STAT/DE User Guide” located at www.globalpointofcare.abbott.
- Caution:** New analyzers or analyzers that have been repaired and returned or replaced will have the factory default settings in the customization profile, as indicated by the DEFAULT0 on the Analyzer Status page. If analyzers in your facility do not use the default customization profile, the appropriate customization profile should be installed before a new, repaired or replaced analyzer is put into use.
- Perform Quality Check** Use the Electronic Simulator to verify the cartridge-reading performance of new or repaired analyzers.

DESCRIPTION

Specifications



| | |
|---------------------------|--|
| DIMENSIONS | Width 7.68 cm (3.035 in.) Length 23.48 cm (9.245 in.) Depth 7.24 cm (2.85 in.) |
| WEIGHT | With rechargeable battery 650 grams (22.9 oz.) With disposable battery 635 grams (22.4 oz.) |
| POWER | Two 9-volt lithium batteries, or rechargeable battery. |
| CALIBRATION | Factory: electronic, mechanical, thermal, pressure |
| MEMORY/CLOCK BACKUP POWER | Lithium Battery |
| DISPLAY | Dot matrix supertwist liquid crystal |
| COMMUNICATION LINK | Infrared light-emitting diode (LED) |
| OPERATING TEMPERATURE | 16-30°C (61-86°F) for i-STAT cartridge testing |
| TRANSPORT TEMPERATURE | -10-46°C (14-115°F) |
| RELATIVE HUMIDITY | 10-90% non-condensing |
| BAROMETRIC PRESSURE | 300-850 mmHg |
| LASER SCANNER | Complies with U.S. 21 CFR 1040.10 and 1040.11 except for deviations pursuant to laser Notice No. 50, dated June 24, 2007. EN 60825-1:2014 IEC 60825-1:2014 |

Software

All analyzer functions are controlled by application software that can be updated as additional tests and features are developed. Coefficients used to maintain the accuracy of cartridge results over time are programmed into the analyzer via CLEW software updates every six months. See under “Standardization and Calibration” in Section 3 of this manual for an explanation of CLEW.

Power

There are two power options for the analyzer: disposable and rechargeable. The analyzer is shipped with a battery carrier for use with two disposable Ultralife 9-Volt lithium batteries (APOC List Number: 06F21-26). Ultralife 9-Volt lithium batteries are manufactured by Ultralife Battery and Energy Products and sold by Abbott Point of Care for use with the i-STAT 1 Analyzer. Only i-STAT rechargeable batteries (APOC List Number: 06F23-55) may be used.

Note: The Ultralife 9-volt lithium battery (APOC List Number: 06F21-26) has a safety feature that provides protection preventing the i-STAT 1 Analyzer from overheating due to component failure within the analyzer circuitry.

Battery Compartment

The battery compartment is located at the display end of the analyzer next to the laser barcode scanner window. The procedure for changing disposable and rechargeable batteries can be found in the Routine Care of the Analyzer and Downloader section of this manual.

Disposable Batteries

The analyzer requires two 9-volt lithium batteries. The lifetime for a set of batteries is mainly dependent on the mix of cartridges in use. Cartridges that require thermal control consume more energy because of heating. Coagulation and immunoassay cartridges consume more energy because of the longer test cycle. A minimum of 400 thermally controlled cartridge uses, about 100 coagulation cartridges, or 50 immunoassay cartridges can be expected before replacement is necessary. Backlighting, if used continuously, may reduce battery life up to 50%. Extensive laser scanning will affect battery life slightly.

The lithium batteries should be removed from the analyzer when long periods, such as six months, of no use are anticipated.

Rechargeable Battery

The analyzer can be powered by a nickel-metal-hydride rechargeable battery. The battery capacity for one full charge is 30% (minimum) of the capacity of one set of disposable lithium batteries (see above). If the analyzer is not in use, batteries will lose approximately 10-30% of their charge over 30 days if not recharged.

Store rechargeable batteries in a cool dry place when not in use.

The battery recharges when the analyzer is placed in a Downloader/Recharger. The battery pack can be removed from the analyzer and placed in the separate recharging compartment on the Downloader/Recharger. Full recharge from a discharged state takes approximately 40 hours. The analyzer will display "Low Battery" when battery recharge is needed.

Caution: Do not short circuit, incinerate or mutilate the rechargeable batteries.

Low Battery Warning


The analyzer will display "Low Battery" when the On/Off key is pressed. Additionally, a flashing battery icon will display on the results screens, as well as the Test Menu and Administration Menu screens when battery replacement is needed. Data is not lost when batteries are fully discharged.

Additional Power


A lithium battery inside the analyzer maintains the clock/calendar and customization profile. This battery should last seven years.

Cartridge Port

Cartridges and the Electronic Simulator are inserted into the analyzer through the cartridge port on the keypad end of the analyzer.

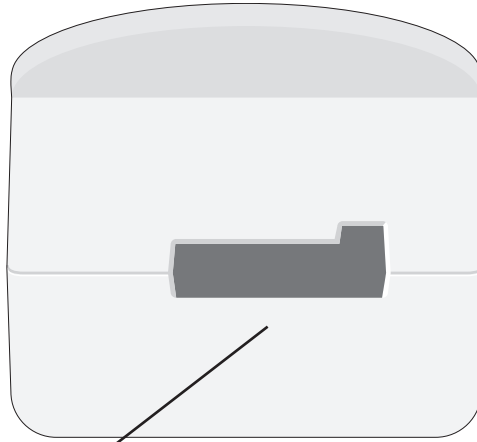

i-STAT Cartridge Port

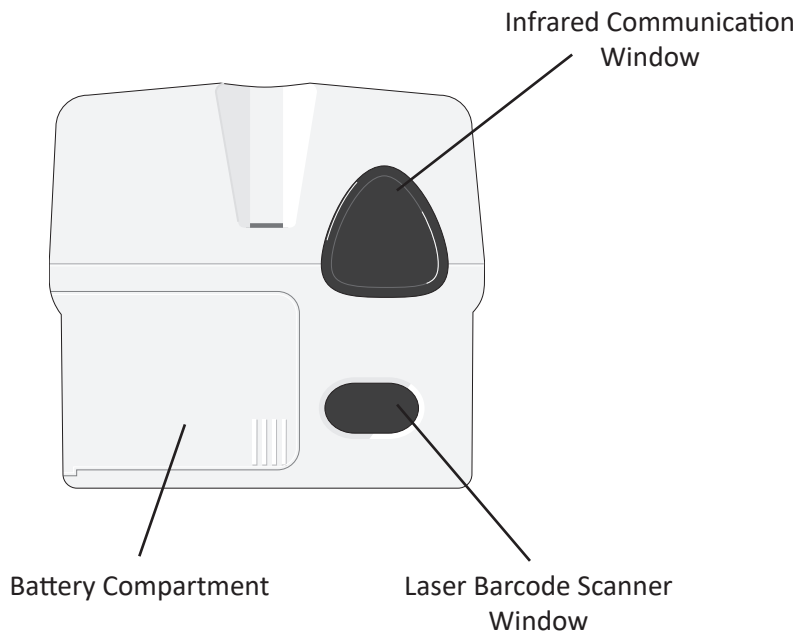

Infrared Communication Window

The Infrared Communication Window provides the analyzer with two-way communication to the Central Data Station via a Downloader, allows analyzer-to-analyzer software updates, and allows analyzer-to-printer communication for printing.

Thermal Control

The analyzer contains a thermal control subsystem of thermistors and heating contact wires that controls the temperature of the sensors and fluids that come into contact with the sensors to 37°C. This subsystem is activated automatically when a cartridge containing tests which require thermal control at 37°C is inserted into the analyzer.

Barometric Pressure Sensor

The analyzer contains a solid-state barometric pressure sensor, which determines the ambient atmospheric pressure used for the PO_2 sensor calibration.

Cartridge Test Cycle

An operator starts a cartridge test cycle either by selecting i-STAT Cartridge from the Test Menu or Quality Tests from the Administration Menu.

The analyzer:

- makes electrical contact with the cartridge
- identifies the cartridge type
- releases calibration fluid to the sensors (when applicable)
- mixes sample and reagent (when applicable)
- measures barometric pressure
- heats the sensors to 37°C (when required by the test)
- measures electrical signals generated by the sensors and calibration fluid (when applicable)
- displaces the calibrant solution with sample (when applicable)
- measures electrical signals generated by the sensors and sample
- accepts the operator and patient IDs scanned or entered by the operator
- accepts chart page information
- calculates and displays results
- stores results

Data Entry

Data that can be scanned into the analyzer or entered via the keypad include:

- Operator ID
- Patient ID, Proficiency ID, or Simulator ID
- Cartridge Lot Number
- Control Lot Number
- Cal Ver Kit Lot Number
- Comment codes for patient and control results
- Chart Page
 - Sample Type
 - Patient Temperature - The analyzer will interpret numbers between 50.0 and 110.0 as degrees Fahrenheit and between 10.0 and 45.0 as degrees centigrade. When a patient temperature is entered, blood gas results will be displayed at both 37°C and the patient's temperature.
 - FIO2
 - Free Fields: three fields, up to 9 characters each

Pt: 145

Scan or Enter Data

Sample Type

Field 1 -----

Field 2 -----

Field 3 -----

PtTemp -----

FIO2 -----

CPB ----- No

1-ART 4-CAP

2-VEN 5-CORD

3-MIX 6-OTHR

→ Page

Results Ready

See the Customization section in this manual for barcode formats recognized by the analyzer.

Storage of Results

The analyzer automatically stores up to 1,000 test records. A test record consists of:

- a set of results
- the date and time the test was performed
- the cartridge type
- all information entered by barcode scanner or keypad including:
 - Operator and Patient IDs
 - Lot numbers for controls and cartridges
 - Chart page data
 - Serial number of the Electronic Simulator
- the serial number of the analyzer
- the number of times the analyzer has been used
- the software and CLEW versions installed in the analyzer
- the name of the analyzer's customization profile

Quality Check Codes, which may appear during the test cycle indicating a problem with the sample, calibration, sensors, mechanical or electrical functions of the analyzer, are also stored.

The Analyzer Status option under the Administration Menu lists the number of stored records as "Total" and "Unsent" records. Test records are stored as "Unsent" until the analyzer uploads data to i-STAT/DE at which time the records are marked as sent. The analyzer can be customized to display a Memory Full prompt or to disable testing until data is transmitted to i-STAT/DE. Otherwise, the oldest data is overwritten when the memory becomes full. Stored test records can be reviewed through the Data Review option on the Administration Menu screen described later in this section.

LCD Display and Backlight


Test results, operator prompts and other messages are displayed on the analyzer's LCD Screen. The backlight for the display is turned on and off by pressing the 0 key for one second. The backlight will automatically turn off after ninety seconds and when the analyzer powers down or is turned off. The backlight cannot be turned on while data entry screens are displayed.

Audible Indicator

The analyzer will beep to indicate:

- whenever a key is pressed.
- a successful barcode entry.
- results are ready.
- a Quality Check Message is displayed.

The analyzer can be customized to disable beeping when a key is pressed or results or messages are displayed.

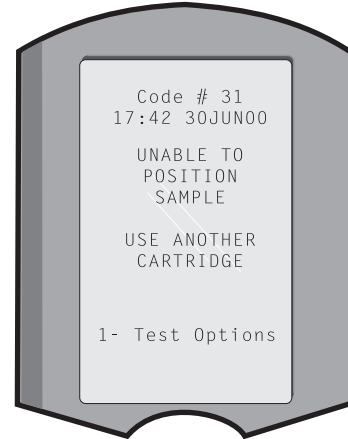

Time Out

The analyzer automatically turns off after a certain period of inactivity.

- **Results displayed:** Results are displayed for 2 minutes before the analyzer turns off provided that a mandatory Comment Code prompt is not displayed. This Inactivity Time Out default time can be increased using Customization.

If a mandatory Comment Code prompt is displayed, the analyzer will turn off after 15 minutes or after the Inactivity Time Out, whichever is greater. In the case of a missed required Comment Code, results will be stored and “_ _ _” will be entered as the Comment Code.










- **Prompting for mandatory data when results are ready for display:** The analyzer will turn off after 15 minutes or after the Inactivity Time Out, whichever is greater, if there is no response to a mandatory data prompt. A mandatory data prompt is a prompt for information that must be entered before pending results are displayed.

In the case of a missed mandatory data prompt, results will not be stored and the test record will state “Test Cancelled by Operator.”

- **Waiting for insertion of cartridge:** After the prompt “Insert Cartridge” is displayed, the analyzer will wait 15 minutes for the operator to insert a cartridge unless the analyzer is in the Proficiency path, in which case the analyzer will wait 5 minutes. If a cartridge is not inserted, the analyzer will turn off. This timeout cannot be customized.
- **Other:** The analyzer will turn off after 2 minutes of inactivity (no keys pressed) in all other circumstances.

Keypad

There are 19 keys located directly below the display. When using the keypad to enter information, the number of dashes in the data entry line will indicate how many characters can be entered on the line. The dash where the next entry will be placed will flash.

| Key | Function |
|---|--|
| SCAN | Activates the barcode scanner. Information that can be entered into the analyzer via the scanner includes: operator ID, patient ID, control and cartridge lot number, patient chart data and comment codes. |
|  | Used to move the cursor on the Set Clock screen and to move up and down the alphabet when the ABC key is pressed. The  (right arrow) key is used as a page key to move from one screen to the next. When Patient ID Recall is enabled, the  key will recall the last patient ID when the analyzer is prompting for Patient ID. The  (left arrow) key is used to backspace and clear keypad entries, and to move backward through the screens within a menu. |
| ABC | Used to enter alpha characters on data entry screens. When the ABC key is pressed the letter A is entered. The arrow keys are used to move up and down the alphabet. To enter a second letter, press the ABC key once to move to the next position and again to enter an A. To enter a number after a letter, press a numbered key. To erase a letter, press the ABC key to move to the next position, then use the  key to backspace and clear the letter. |
| 0 – 9 | Used to enter digits on data entry screens and to select menu options and stored records. |
| . | Enters a decimal point or a comma separator according to the analyzer's Customization Profile. |
|  | Used to turn the screen backlight on and off. |
| Enter  | Used to respond to a prompt to complete an action, such as entering an operator or patient ID via the keypad. |
| MENU | Used to return to the previous menu and switch between the Test and Administration Menus. |
| Print  | Used to print either directly to the portable printer or to the portable printer attached to a Downloader. |
| On/Off  | Turns the analyzer on or off. When the analyzer is on, the On/Off key must be pressed for a second to turn the analyzer off. This key is inactive when a test is in progress and when the analyzer is prompting for mandatory data. |

i-STAT 1 Menu Tree There are two main menus: The Test Menu and the Administration Menu.

| Test Menu | Administration Menu | |
|---------------------|---------------------|--|
| 1- Last Result | 1. Analyzer Status | Temp Pressure Battery Uses Serial CLEW Version Custom Stored Records Total Unsent |
| 2- i-STAT Cartridge | 2- Data Review | 1-Patient 2-Control 3-Proficiency 4-Cal Ver 5- Simulator 6- All 7- List |
| | 3-Quality Tests | 1-Control 2- Proficiency 3- Cal Ver 4- Simulator |
| | 4- Customization | 1-View 1- Analyzer 2- ID Entry 3- Patient Tests 4- QC Tests 5- Results 2-Change 1- Analyzer 2- ID Entry 3- Patient Tests 4- QC Tests 5- Results 6- Password 7- Restore Factory Settings |
| | 5- Set Clock | |
| | 6- Transmit Data | 1- Most Recent 2- This Month 3- Last Month 4- All 5- Unsent |
| | 7-Utility | 1- Send Software 2- Clear Memory 3- Receive Software |

TEST MENU



The Test Menu is displayed when the analyzer is turned on using the On/Off key.

The options are:

- 1 - Last Result
- 2 - i-STAT Cartridge

Option 2 is used for testing patient samples.

Note: If the handheld is customized to disable testing under certain conditions, the disabled option will be listed without its number so that it cannot be selected.

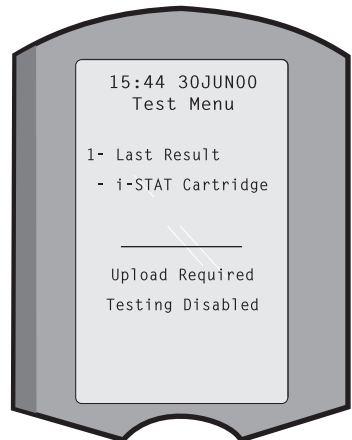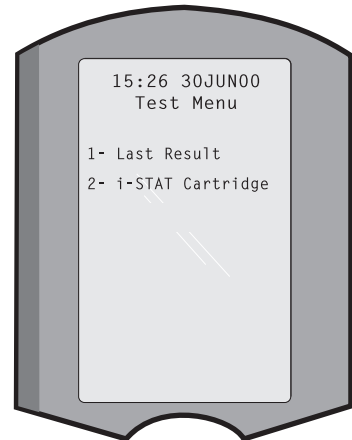


ADMINISTRATION MENU

Overview

The Administration Menu is accessed by pressing the Menu key from the Test Menu screen. The options are:


- 1 - Analyzer Status
- 2 - Data Review
- 3 - Quality Tests
- 4 - Customization
- 5 - Set Clock
- 6 - Transmit Data
- 7 - Utility

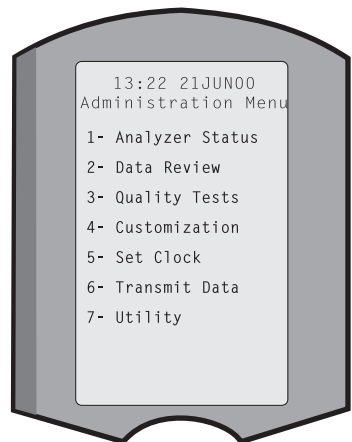

Analyzer Status

The Analyzer Status screen contains information about the condition or “status” of the analyzer. Fresh readings are made whenever this option is selected.

| | |
|-----------------|---|
| Temp | Room temperature. |
| Pressure | Barometric pressure. |
| Battery | Battery voltage. |
| Uses | Total number of cartridge and simulator test cycles, whether or not results reported. |
| Serial | Serial number of the analyzer. |
| CLEW | Version of standardization data installed in the analyzer. |
| Version | Version of application software installed in the analyzer. |
| Custom | Customization profile name. |

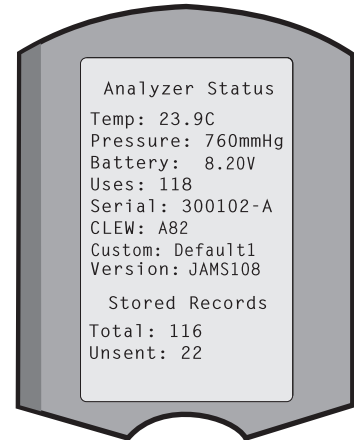

| | |
|-----------------------|---|
| Stored Records | Total: The number of test records in the analyzer’s memory. The maximum storage capacity is 1,000 test records, which include records with results and Quality Check Codes for patients and controls both liquid and electronic. Unsent: The number of test records that have not been transmitted to i-STAT/DE. |
|-----------------------|---|

Data Review

The Data Review function allows the operator to review stored results by the categories listed below. The number of test records stored is indicated at the bottom center of the screen as x/y where x is the record on the screen and y is the total number of stored records in the selected category. The 1 and 2 keys are used to scroll through the stored records as indicated on the bottom right and left of the screen. The most recent test record is always in the first position. The right arrow key is used to page through the screens of the displayed record.

1 - Patient The records for a patient are recalled by scanning or entering via the keypad the Patient ID. If no Patient ID is entered, all patient tests are recalled.


2 - Control

3 - Proficiency

4 - Cal Ver


5 - Simulator All external and internal Electronic Simulator records.

6 - All All test records in the analyzer’s memory.

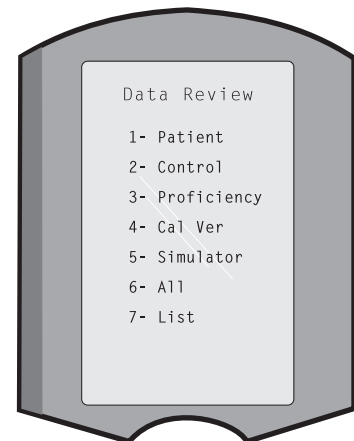

7 - List

Records are listed with Cartridge type, date and time of test, patient ID, control lot, proficiency ID, or Cal Ver lot and test level as applicable. Any number of test records can be selected for viewing or printing using the number keys. Pressing the number key corresponding to a record selects a record; pressing the number key a second time deselects the record.

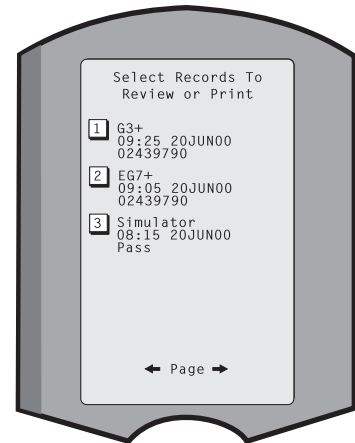

To view one or more records, select the records and press the Enter key. To print records, select the records and press the Print key.

Quality Tests

Non patient tests can be initiated from the Quality Tests menu. Options are:

- 1 - Control
- 2 - Proficiency (external quality control)
- 3 - Cal Ver (Calibration Verification for cartridges)
- 4 - Simulator (cartridge-reading function only)

When testing is initiated from one of these options, the handheld prompts the operator to scan or enter the Operator ID; the Control Lot Number, Proficiency ID, Cal Ver Kit Lot Number, or Simulator ID as applicable; and the Cartridge Lot Number.

When the Quality Tests option is used, results can be reviewed according to the corresponding options under the Data Review option.



Customization

Analyzers can be customized for site-specific testing characteristics and requirements. A complete list of customizable parameters and their default values can be found in the Customization section. An analyzer can be customized via the keypad or via i-STAT/DE. Items that cannot be customized via the analyzer's keypad are operator lists, test strip lists, reference and action ranges, sample types and order of items on the Chart page.

i-STAT/DE's Customization function can be used to create one customization profile for all analyzers or different profiles for different locations. When the Customization function is enabled, the profiles are transmitted to the analyzers when they are placed in a Downloader.

Caution: If location specific customization profiles are created, analyzers should not be moved from one location to another unless they are re-customized for the new location. This is especially important if "CPB: Automatically Adjust" or "CPB: Do Not Adjust" is included in a location-based customization profile. The CPB function adjusts hematocrit and hemoglobin results for the dilutional affect of pump fluid during cardiopulmonary bypass surgery. If an analyzer customized for the CVOR as "CPB: Automatically Adjust" is used for patients who are not on the pump, hematocrit results will be reported falsely high. If an analyzer customized as "CPB: Do Not Adjust" is used for patients who are on the pump, hematocrit results will be reported falsely low. For details on the CPB function, see the Theory section of this manual.

It is recommended that only one method, i-STAT/DE or the keypad, be used to customize all analyzers within a site. If both methods are in use, and the Customization function is not disabled in i-STAT/DE, any changes made to the profile of an analyzer via the keypad will be overwritten the next time the analyzer is placed in the downloader.

The customization profile of an analyzer is identified in the Customization option under the Administration Menu on the analyzer. DEFAULT0 indicates that the analyzer has factory settings. When an analyzer has been customized via i-STAT/DE, the name assigned to the profile by i-STAT/DE is listed. If the default or i-STAT/DE profile is changed on the analyzer, the profile is listed as 00000000.

Viewing the Customization Profile

Select **4- Customization** from the Administration Menu, select **1- View** then select from the Customization Menu:

- 1 - Analyzer
- 2 - ID Entry
- 3 - Patient Tests
- 4 - QC Tests
- 5 - Results

Select a category to review. Use the ← and → keys to scroll through the preferences for each category and use the ← key to return to the Customization menu.

The Customization review option on the analyzer does not display the certified operator list. This item can be viewed in i-STAT/DE.

Note:

- Outside the USA, the following changes should be considered: language, unit set, date format and decimal separator.

1 - Analyzer

First page

Date Format
Sound
Auto-transmit
Memory Full
Batch Mode Timeout

Second page

Inactivity Timeout
Upload Schedule
Clock Password
Sync Clock
Patient Record Limit Access

Third page

Wireless (only available with the i-STAT wireless analyzer)

2 - ID Entry

1 – Operator ID

First page

Minimum Length
Maximum Length
Repeat ID
Manual Entry
Code I2of5

Second page

Code 128
EAN-8, EAN-13
Codabar
Code 93
Code 39

Third page
Code 39 Check Digit
Truncate First
Truncate Last
Operator List
Not Certified Action
Not In List Action

Fourth page
Warn User
Print ID

2 – Patient ID

First page
Minimum Length
Maximum Length
Repeat ID
ID Recall
Manual Entry

Second page
Code I2of5
Code 128
EAN-8, EAN-13
Codabar
Code 93

Third page
Code 39
Code 39 Check Digit
Truncate First
Truncate Last
Patient List
Not in List Action

Fourth page
Lockout Override
Confirmation Method
Print ID

3 - Patient Tests

First page
Cartridge Auto-chart
Cartridge Information
Cartridge Barcode
Cartridge Lot Number
Comment Code In Range

Second page
Comment Code Out of Range
Cart Sample Type
Result Output
Downloader Lockout
STATNotes

- 4 - QC Tests
 - 1 – Simulator
 - Ext Simulator
 - Int Simulator
 - Int Simulator Schedule Option
 - 2 – Cartridge QC
 - [First page](#)
 - Pass/Fail Method
 - Comment Code In Range
 - Comment Code Out of Range
 - Result Format
 - APOC fluid Lot Scan Only
 - [Second page](#)
 - eVAS Name
- 5 - Results
 - 1 – ACT/Ref Ranges
 - 2 – Display Ranges
 - 3 – Units
 - 4 – Options
 - [First page](#)
 - Decimal Separator
 - Test Selection
 - Hematocrit
 - Base Excess
 - ACT-C
 - [Second page](#)
 - ACT-K
 - Print Ref. Ranges

Changing the Profile

To customize via the handheld keypad, select **4- Customization** from the Administration Menu, then select **2- Change**. If the handheld has already been customized with a password, enter the password. If not, press the Enter key. (It is recommended that the Change function be password protected). Then make selections from the Customization menu. To change a setting, select the item by pressing the number key corresponding to the item, then select the setting. Use the → key to view all items. After all items have been set, turn the handheld off to save and activate the settings.

Note:

- Outside the USA, the following changes should be considered: language, unit set, date format and decimal separator.

- 1 - Analyzer
 - [First page](#)
 - Language
 - Date Format
 - Sound
 - Auto-transmit
 - Memory Full
 - [Second page](#)
 - Batch Mode Timeout
 - Inactivity Timeout
 - Upload Schedule
 - Clock Password
 - Sync Clock
 - [Third page](#)
 - Wireless (available with the i-STAT 1 wireless)

2 - ID Entry

1 – Operator ID

First page

Minimum Length
Maximum Length
Repeat ID
Manual Entry
Code I2of5

Second page

Code 128
EAN-8, EAN-13
Codabar
Code 93
Code 39

Third page

Code 39, Check Digit
Truncate First
Truncate Last
Print ID

2 – Patient ID

First page

Minimum Length
Maximum Length
Repeat ID
ID Recall
Manual Entry

Second page

Code I2of5
Code 128
EAN-8, EAN-13
Codabar
Code 93

Third page

Code 39
Code 39 Check Digit
Truncate First
Truncate Last

3 - Patient Tests

First page

Cartridge Auto-chart
Cartridge Information (functionality preset by analyzer firmware)
Cartridge Barcode (functionality preset by analyzer firmware)
Cartridge Lot Number (functionality preset by analyzer firmware)
Comment Code, In Range

Second page

Comment Code, Out of Range
Cart Sample Type
Result Output
Downloader Lockout

4 - QC Tests

1 – Simulator

Ext Simulator
Int Simulator
Int Simulator Schedule Option

- 2 – Cartridge QC
 - Pass/Fail Method
 - Comment Code In Range
 - Comment Code Out of Range
 - Result Format
 - APOC Fluid Lot Scan Only

- 5 - Results

- 1 – Units and Ranges

- 2 – Options

- [First page](#)

- Decimal Separator

- Test Selection

- Hematocrit

- Base Excess

- ACT-C

- [Second page](#)

- ACT-K

- Print Ref. Ranges

- 6 - Password

- 7 - Restore Factory Settings


Note: For additional procedures related to customization using the analyzer keypad, refer to *Procedures for Customization using the Analyzer Keypad* in Section 9 of this manual.

Set Clock

If the analyzer is customized with a password, the Set Clock function will be password protected. If a password has not been assigned, pressing the Enter key will display the time and date screen. Use the arrow keys to move the cursor to the digit to be changed. Use a number key to change the digit. Press Enter to accept the changes or Menu to cancel the changes. An invalid entry, such as 13 for a month, will not be accepted.

The format of the date on this screen can be customized using the i-STAT/DE customization function, as mm/dd/yy or dd/mm/yy. The analyzer recognizes years in which February has 29 days.

The analyzer can be customized using i-STAT/DE to synchronize or update the real time clock to the i-STAT/DE's clock at the time of each download. This option eliminates the need to reset the analyzer's clock at the beginning and end of Daylight Saving Time. Otherwise, the clock must be manually changed for Daylight Savings Time changes.

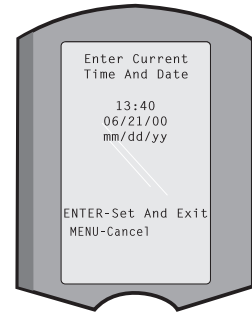

Transmit Data


Unsent test records are automatically transmitted to i-STAT/DE when an analyzer is placed in a Downloader/Downloader/Recharger. In some cases it may be desirable to have the capability to retransmit data. The Transmit Data function allows transmission of data in the following manner:

- 1 – Most Recent
- 2 – This Month
- 3 – Last Month
- 4 – All
- 5 – Unsent

Most Recent is the result from the last cartridge tested.

The analyzer can be customized using i-STAT/DE to apply a date range limit to the Transmit All functions.


Auto-transmit is temporarily disabled when the Transmit Data option is selected to allow the user to control transmission of data.

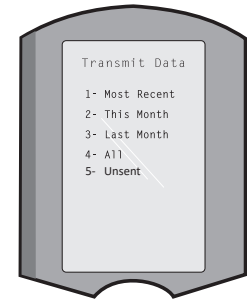

Utility

The Utility menu can be password protected using the Customization function on the analyzer or in i-STAT/DE.

- 1 – Send Software: Allows the analyzer to transmit software to another analyzer. See the Software Update section of this manual.
- 2 – Clear Memory: Erases results from the analyzer's memory. Options are:
 - 1 – Previous to 01MMYY (where MMYY is current month and year, such as 01JUN00)
 - 2 – Previous to 01mmmyy (where mmmyy is previous month and year, such as 01May00)
 - 3 – All
 - 4 – Cancel
- 3 – Receive Software: Allows users to remotely request a JAMS and CLEW update for the analyzer from i-STAT/DE. See section 18 (Updating Software) for full details.

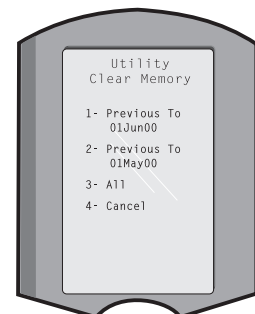

Caution

Do not open the analyzer. The analyzer may only be opened by factory authorized service personnel. **Class 2 laser radiation when open; DO NOT stare into the laser aperture or the laser beam, or point the laser beam at other persons.**

Use of controls, adjustments or performance of procedures other than those specified herein may result in hazardous laser radiation exposure.

Class 2 laser scanners use a low power, visible light diode. As with any bright light source, such as the sun, the user should avoid staring directly into the laser beam. Momentary exposure to a Class 2 laser is not known to be harmful.

Barcode Label Quality

To ensure that printed barcode labels are reliably read by i-STAT handhelds, the best available printing methods and settings should be used. However, as specified in the *Health Industry Bar Code (HIBC) Provider Applications Standard* (ANSI/HIBC 1.3-2010), the quality of printed labels should meet the minimum grade level of 1.5.

Ambient Lighting from LED Light Sources

The analyzer's barcode scanning functionality may experience interference when scanning barcodes under ambient light from an LED light source. This interference results in the analyzer being unable to scan a barcode at all (no beep acknowledgement). When scanning barcodes under ambient light from an LED light source, it is recommended that the barcode be shielded from the ambient light when attempting to scan the barcode.

Procedure

Before scanning, check to see what information is required by the displayed prompt. Hold the analyzer 3-9 inches (8 – 23 cm) from the barcode to be scanned. An angle of about 10 degrees from perpendicular is best. Hold the analyzer and place the object to be scanned on a flat surface or, place the analyzer on a flat surface and hold the object in front of the analyzer. Avoid accidentally scanning other nearby items. Avoid pointing the beam into anyone's eyes.

STEP

ACTION

-
- | | |
|---|---|
| 1 | Press and hold down the Scan key to start the barcode scanner. The analyzer emits a visible red beam. |
| 2 | Position the analyzer and barcode so the beam forms a red line that spans the entire barcode. Increasing distance between the barcode and analyzer lengthens the red line. The analyzer does not need to touch the barcode. |
| 3 | When the analyzer accepts the barcode, it will beep in acknowledgement and automatically turn off the beam. The beam will also turn off after 3-4 seconds. |
| 4 | View the data that was scanned by the analyzer and verify that it is correct. |
| 5 | Release the Scan key. |

Note: If the Scan key is released as soon as the beep is heard, the next prompt will be displayed and the information scanned will not be able to be viewed.

PROMPTS AND MESSAGES

Prompts Either before or during the testing cycle, the analyzer will display prompts that require an operator action or keypad entry, such as “Enter Operator ID.” Prompts are described in the manual when used. Some prompts require input before results are displayed. Prompts for the following information are mandatory:

- Operator ID
- Patient ID
- Lot Numbers for Quality Tests
- Cartridge Lot Number

Startup Messages When the On/Off key is pressed the analyzer may display one or more startup messages. A startup warning message indicates an action that should be taken in the near future to maintain the analyzer in working condition. If the analyzer is customized to disable testing under certain conditions, a startup lockout message indicates the action that must be taken before testing is re-enabled.

Quality Check Messages If the analyzer detects a problem during power on, a Quality Check message will be displayed indicating the action that must be taken before testing can begin.

A Quality Check message will also be displayed and testing halted if the analyzer detects a problem during the test cycle.

Startup messages and Quality Check messages are described in the Troubleshooting section of this manual. “Upload Required, Testing Disabled” is an example of a startup lockout message, “Battery Low” is an example of a startup warning message, and “Unable to Position Sample” is an example of a quality check failure during the testing cycle.

Note: The “Cartridge Locked” or “Simulator Locked” prompt is always displayed when a cartridge or Electronic Simulator is inserted into the analyzer. Any attempt to remove a cartridge or Electronic Simulator before this prompt is removed from the screen may cause damage to the analyzer.

© 2022 Abbott. All rights reserved. All trademarks referenced are trademarks of either the Abbott group of companies or their respective owners.



i-STAT 1 SYSTEM TECHNICAL BULLETIN

NETWORK OPTIONS FOR UPDATING THE i-STAT 1 ANALYZER USING www.globalpointofcare.abbott

OVERVIEW

This Technical Bulletin has been created specifically to guide you through the process of updating the software on the i-STAT 1 Analyzer(s) using one of three network options:

1. Using a Network Downloader, Downloader/Recharger, or DRC-300 and the JammLite process with TCP/IP
2. Using a Serial Downloader or a serially connected DRC-300 and the JammLite process for accounts with RALS-Plus and i-STAT/DE
3. Using the i-STAT/DE Customization Workspace and a Network Downloader, Downloader/Recharger, or DRC-300

Before beginning, check the Analyzer Status page and verify that the i-STAT 1 Analyzer has enough battery power (7.5 volts or higher).

WHY DO i-STAT ANALYZERS NEED SOFTWARE UPDATES?

The i-STAT System is designed to eliminate operator influence on delivered results.

Due to continuous manufacturing process improvements to the i-STAT System, it is necessary to update standardization values from time to time to maintain long-term consistency of performance. These updates are equivalent to manually adjusting calibration on a traditional laboratory analyzer.

New CLEW software -- delivered twice a year -- re-establishes these standardization values and incorporates refinements to the internal quality monitoring system. New JAMS application software allows the i-STAT Analyzer to recognize any newly launched cartridge types and to perform any newly launched features.

i-STAT is for *in vitro* diagnostic use.

1. UPDATE PROCEDURE USING A NETWORK DOWNLOADER AND THE JammLite PROCESS WITH TCP/IP

1.1: Before starting the process, make sure all the required equipment / information is available.

• **Computer with:**

- Windows 2000, XP, or Windows 7
- Access to www.globalpointofcare.abbott

• **i-STAT System Equipment**

- (1) Network Downloader, Downloader / Recharger, or DRC-300

Note: These instructions assume that the Network Downloader types used for the update process are already installed and in use on the customer's network.

- (2) i-STAT 1 Analyzer
- (3) Electronic Simulator

• **List of the IP Address(es) for the network downloader(s) to be used for the software update process**

1.2: Close all open programs on the computer.

1.3: Navigate to www.globalpointofcare.abbott > Support > i-STAT 1 and i-STAT Alinity Support > i-STAT 1 Resources Login > Product Software > i-STAT System Software Updates and Access Software.

1.4: Scroll to "Step 2: Download Software Update File".

No Data Management

Click on "DOWNLOAD SUXXXXXX.ZIP" and save the file to the Desktop. Close the "Download Complete" window.


Data management with i-STAT/DE version ≥2.8.0.1

Click on "DOWNLOAD SUXXXXXX.ZIP" and save the file to the Desktop. Close the "Download Complete" window.

Note: i-STAT/DE versions less than 2.8.0.1 are no longer supported. Use a Network Downloader, Downloader/Recharger, or DRC-300 and the JammLite process to perform the software update via its IP Address (TCP/IP) or Port (COM Port) for serially connected downloaders.


1.5: Navigate to saved zip file location. Right click on the zip file and select Extract All and Extract to the Desktop.

- Navigate to the Desktop and click on the folder SUXXXXXX to open.


1.6: Double click the software file “SUXXXXXX.exe” to run. If a Command window opens prompting to overwrite, answer “Y” and then press Enter. Continue answering “Y” to all prompts that appear until the Command window closes. From among the icons that appear, double click  JAMMLITE.EXE to launch the JammLite Utility.

- If the JammLite program does not launch or you receive an error message, contact APOC Technical Support and tell the support specialist you are unable to launch the JammLite Utility.


1.7: In the JammLite utility, select the **i-STAT 300 Analyzer** within the Instrument dropdown menu.

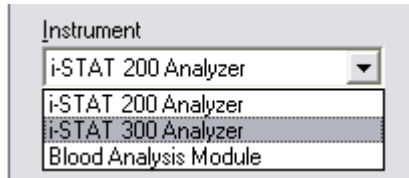

1.8: Select **TCP/IP** within the Port dropdown menu.

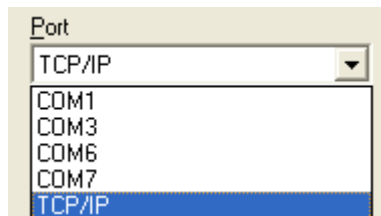


1.9: Type the IP Address of the Network Downloader being used for the software update in the **IP Address** box.

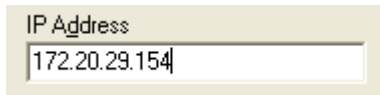


Note: the address used is for example only.

1.10: Check that the **Application** and **CLEW** listings match those in the Product Update. Click the **Update** button.

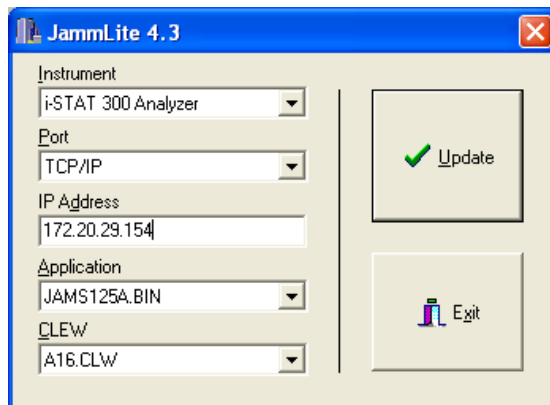


Note: Application and CLEW numbers are for example only.

1.11: Follow the onscreen instructions.

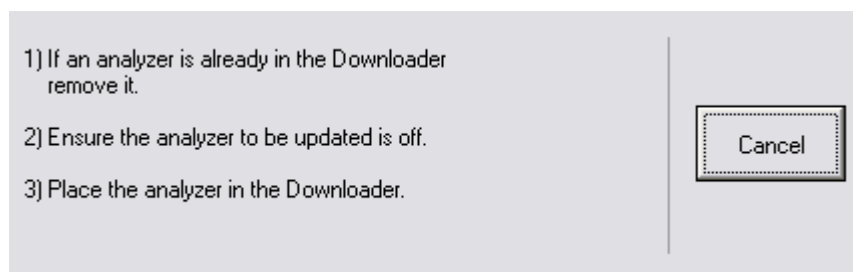


1.12: When the update is in progress, the following screen will appear:

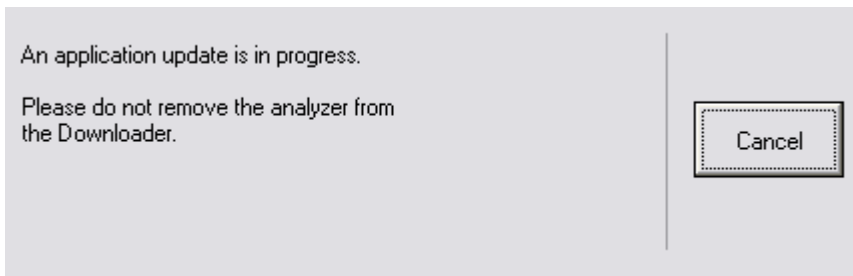

The analyzer will have 1's and 0's streaming across the screen signifying that it is receiving the software.

Do not move the analyzer until the success screen is displayed.

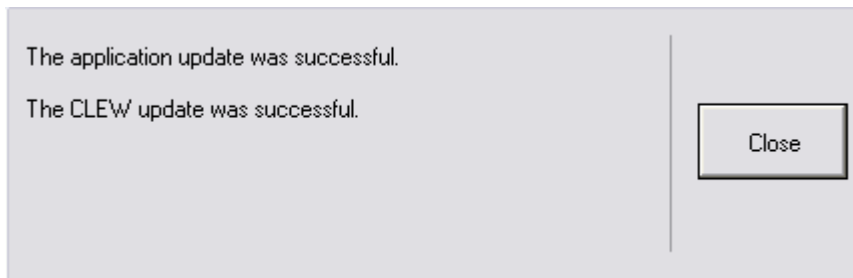


1.13: Run the Electronic Simulator in the analyzer. When the simulator finishes, PASS should be displayed.

Note: If PASS is not displayed, re-run the Electronic Simulator. If the repeated Electronic Simulator attempt fails, contact APOC Technical Support. For additional information on running the Electronic Simulator, see Section 14 of the i-STAT 1 System Manual.

Congratulations. The process for updating the first i-STAT 1 Analyzer is complete.

Review the options below for additional instructions.

If there are no additional analyzers to update, the process is complete.

- Click the  button in the upper right corner of the software screen.
- Close all other open boxes.
- Confirm all messages.

If there are additional analyzers to update via the same Network Downloader address:

- Click **Close**.
- Repeat steps **1.10** through **1.13**

If there are additional analyzers to update via a different Network Downloader address:

- Click **Close**.
- Repeat steps **1.9** through **1.13**

2. UPDATING THE i-STAT 1 ANALYZER USING A SERIAL DOWNLOADER OR SERIALLY CONNECTED DRC-300 AND THE JammLite PROCESS FOR ACCOUNTS WITH RALS-Plus AND i-STAT/DE

2.1: Before starting the process, make sure all the necessary required equipment is available.

- Computer with:
 - RALS-Plus IMS
 - Applicable ports, depending on the type of downloader being used
 - Access to www.globalpointofcare.abbott

- i-STAT System Equipment

- (1) i-STAT 1 Analyzer
- (2) Serial Downloader, Downloader / Recharger, or DRC-300

Note: These instructions assume that the Serial Downloader types being used for the update process are already installed and in use at the customer's facility.

- (3) All supplied downloader cables
- (4) Electronic Simulator

2.2: Using the supplied downloader cables, connect the i-STAT 1 serial Downloader, Downloader/ Recharger, or DRC-300 to an accessible port on the back of the computer, and connect the power supply to a wall outlet or power strip.


2.3: Close all open programs on the computer.

2.4: Disable the RALS remote connections.


- Right click on the red RRC dot in the lower right corner of the display next to the clock and go to "Administrator Settings."
- Enter the RRC Password of the Day. Call 877-627-7257 to obtain the password, or go to <https://www.rals.com/us/home/rals-system/software/rals-password.html> and type "richmond" as the key code. Note: Abbott Point of Care Inc. recommends changing the default password.



- Double click on "COM1."

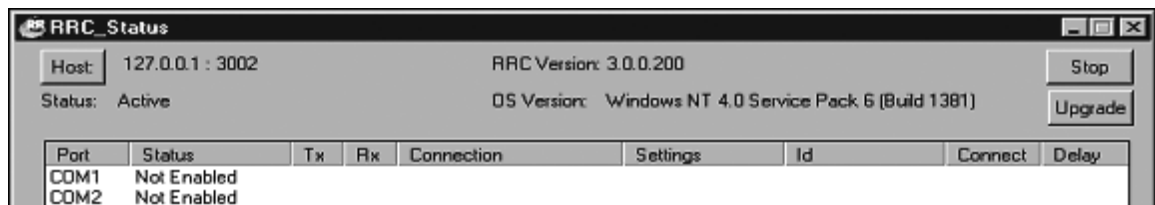

- Record "Host" and "ID" entry.
- Uncheck the "Enabled" box. Click OK.

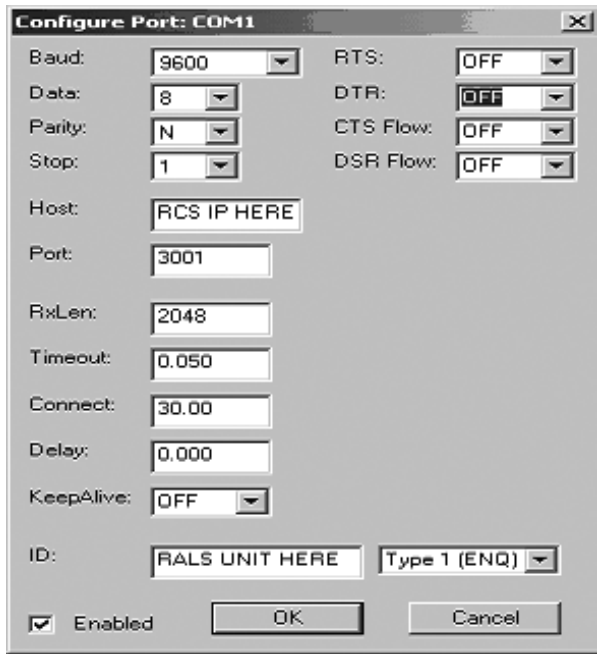

- If not enabled, minimize the window and proceed to the next step.


2.5: Navigate to www.globalpointofcare.abbott > Support > i-STAT 1 and i-STAT Alinity Support > i-STAT 1 Resources Login > Product Software > i-STAT System Software Updates and Access Software.

Scroll to “Step 2: Download Software Update File”




Note: i-STAT/DE versions less than 2.8.0.1 are no longer supported. Use a Network Downloader, Downloader/Recharger, or DRC-300 and the JammLite process to perform the software update via its IP Address (TCP/IP) or Port (COM Port) for serially connected downloaders.

2.6: Navigate to saved zip file location. Right click on the zip file and select Extract All and Extract to the Desktop.


- Navigate to the Desktop and click on the folder SUXXXXXX.exe to open.
- Click Run. Double click the software file “SUXXXXXX.exe” to run. If a Command window opens prompting to overwrite, answer “Y” and then press Enter. Continue answering “Y” to all prompts that appear until the Command window closes.
- Double click  JAMMLITE.EXE to launch the JammLite Utility.

Note: If the JammLite program does not launch or you receive an error message, contact APOC Technical Support and tell the support specialist you are unable to complete step 2.6 of this document.

2.7: In the JammLite utility, select the i-STAT 300 Analyzer within the Instrument dropdown menu.

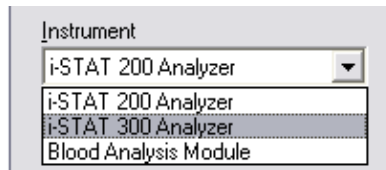


2.8: Make sure there is at least one port number listed under the Port Drop-Down List.

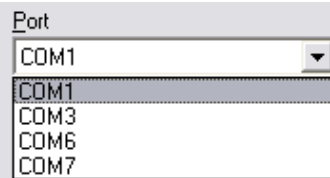

Note: If the Port drop-down list says “None.” make certain the connection from the Serial Downloader to the computer is tight. Then reboot the computer and return to step **2.3**.


2.9: Check that the Application and CLEW listings match those in the Product Update for the current software release.

2.10: Click on the “Update” button.

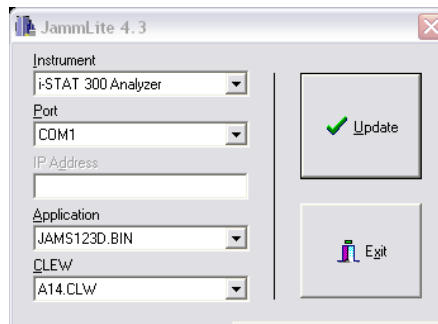


Note: Application and CLEW numbers are for example only.

2.11: Follow the onscreen instructions.

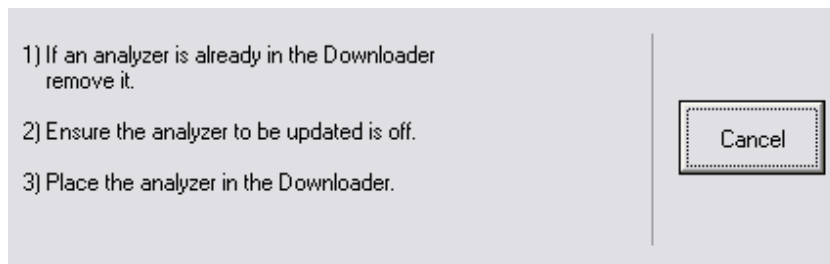


2.12: When the update is in progress, the following screen will appear:

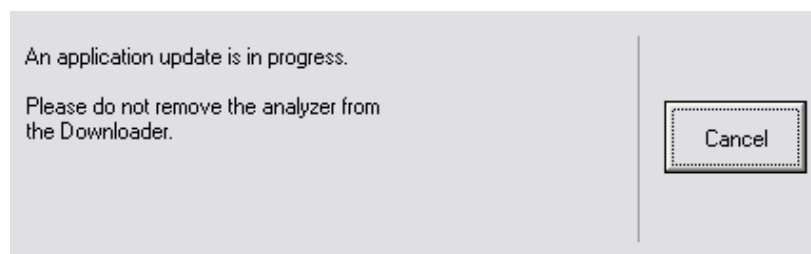

The analyzer will have 1's and 0's streaming across the screen signifying that it is receiving the software.

Do not move the analyzer until the success screen is displayed.

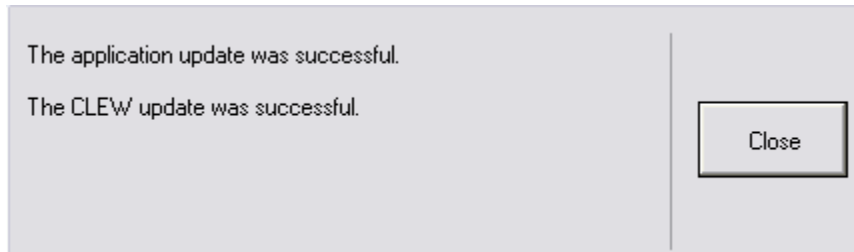

2.13: Run the Electronic Simulator in the analyzer. When the simulator finishes, **PASS** should be displayed.


Note: If **PASS** is not displayed, re-run the Electronic Simulator. If the repeated Electronic Simulator attempt fails, contact APOC Technical Support. For additional information on running the Electronic Simulator, see Section 14 of the i-STAT 1 System Manual.

Congratulations. The process for updating the first i-STAT 1 Analyzer is complete.

If there are additional analyzers to be updated via JammLite:

- Click **Close**.
- Repeat steps **2.10 through 2.13**.



If there are no additional analyzers to update, the process is complete.

- Click the  button in the upper right corner of the software screen.
- Proceed to step **2.14** to restart the RRC if it was disabled in Step **2.4**.
- Proceed to step **2.14** to update the RALS-Plus i-STAT/DE with the latest CLEW and JAMS software files.

2.14: Restart the RALS remote connection.

- Maximize the RRC Status window.
- Double click on Com1.
- Check the “Enabled” box.
- Enter the recorded “Host” and “ID” entry from step **2.4**.
- Click “OK.”

2.15: Update the CLEW version in the Customization Workspace.

- Transfer the files.
 - Access the main Customization Workspace page.
 - Click **Update i-STAT/DE** → **Upload Update File**.
 - Browse to desktop, click on  **XXX.CLW** and click **Upload**.
(Note: the XXX is the CLEW version that you are updating.)
 - Click **Update i-STAT/DE** → **Upload Update File**.
 - Browse to Desktop, click on  **JAMSXXXX.BIN** and click **Upload**.
(Note: the XXXX is the JAMS version that you are updating.)
- In the Customization Workspace, under the “Default customization profile:” column, click on the “**i-STAT Analyzer CLEW**” button.

Note: Customization screens may vary depending upon the i-STAT/DE version in use.

Default customization profile:

Language: English

Unit Set: UNITSET00

i-STAT Analyzer CLEW: A16

Philips BAM CLEW: [None]

i-STAT 1 Software: JAMS125A.BIN

Preferences: DEFAULT0

STATNotes: CHART0

Use Operator List

- Check the box next to the new version of **CLEW** and click **OK**.

Institution: Inst1

Location: Default customization profile

i-STAT Analyzer CLEW:

| Selected CLEW | Name | Expiration Date |
|-------------------------------------|------|----------------------|
| <input checked="" type="checkbox"/> | A16 | 6/24/2009 8:00:00 AM |

Answer OK to the question that appears.

- If “Uses Default” is not checked beside any Location-based customization profile, click the box under the “i-STAT Analyzer CLEW” column.

Location-based customization profiles:

| Location | Enabled | Uses Default | Update CLEW | i-STAT Analyzer CLEW | Philips BAM CLEW | Preferences | STATNotes | <input type="checkbox"/> |
|----------|--------------------------|-------------------------------------|-------------------------------------|----------------------|------------------|-------------|-----------|--------------------------|
| ER | <input type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | A16 | [None] | DEFAULT0 | CHART0 | <input type="checkbox"/> |
| ICU | <input type="checkbox"/> | <input type="checkbox"/> | <input checked="" type="checkbox"/> | A16 | [None] | DEFAULT0 | CHART0 | <input type="checkbox"/> |
| Lab | <input type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | A16 | [None] | DEFAULT0 | CHART0 | <input type="checkbox"/> |
| OR | <input type="checkbox"/> | <input type="checkbox"/> | <input checked="" type="checkbox"/> | A16 | [None] | DEFAULT0 | CHART0 | <input type="checkbox"/> |

- Click the new version of **CLEW**, and then click **OK**.

Institution: Inst1

Location: Default customization profile

i-STAT Analyzer CLEW:

| Selected CLEW | Name | Expiration Date |
|-------------------------------------|------|----------------------|
| <input checked="" type="checkbox"/> | A16 | 6/24/2009 8:00:00 AM |

Answer OK to the question that appears.

- Under the “Default Customization profile:” column, click on the **i-STAT 1 Software** drop-down list. Select the JAMS version that matches the Product Update and click **OK**.

Default customization profile:

Language:
English

Unit Set:
UNITSET00

i-STAT Analyzer CLEW:
A16

Philips BAM CLEW:
[None]

i-STAT 1 Software:
JAMS125A.BIN

Preferences:
DEFAULT0

STATNotes:
CHART0

Use Operator List

3. UPDATING THE i-STAT 1 ANALYZER USING THE i-STAT/DE CUSTOMIZATION WORKSPACE AND A NETWORK DOWNLOADER, DOWNLOADER/RECHARGER, or DRC-300

3.1: Before starting the process, make sure all the required equipment / information is available.

- Computer which can access the Customization Workspace
- Access to www.globalpointofcare.abbott
- i-STAT System Equipment
 - (1) i-STAT 1 Analyzer
 - (2) Network Downloader, Downloader/Recharger, or DRC-300

Note: These instructions assume that the Network Downloader types being used for the update process are already installed and in use on the customer’s network.

- (3) Electronic Simulator

3.2: Update the CLEW and JAMS versions in the Customization Workspace.

- a) Navigate to: www.globalpointofcare.abbott > Support > i-STAT 1 and i-STAT Alinity Support > i-STAT 1 Resources Login > Product Software > i-STAT System Software Updates and Access Software.
- b) Scroll to “Step 2: Download Software Update File”

**Data management with
i-STAT/DE version ≥2.8.0.1**

Click on “**DOWNLOAD SUXXXXXX.ZIP**” and save the file to the Desktop. Close the “Download Complete” window.

Note: i-STAT/DE versions less than 2.8.0.1 are no longer supported. Use a Network Downloader, Downloader/Recharger, or DRC-300 and the JammLite process to perform the software update via its IP Address (TCP/IP) or Port (COM Port) for serially connected downloaders.

- c. Navigate to saved zip file location. Right click on the zip file and select Extract All and Extract to the Desktop.
- d. Access the main Customization Workspace page.
- e. Click **Update i-STAT/DE** → **Upload Update File**
- f. Browse to the Desktop, click on the SUXXXXXX folder to open. Select SUXXXXXX.exe and click upload.

(Note: the XXXXXX is the JAMS and CLEW version that you are updating.) Immediately after uploading the CLEW/JAMS to the i-STAT/DE server, i-STAT/DE will unpack the files and make them available for use in the analyzer’s Customization Workspace.



3.3: Close Windows Explorer by clicking on the  in the upper right corner.


3.4: Access the Customization Workspace.

- RALS-Plus Users:
 - From the RALS-Plus Application, pick i-STAT from the drop-down menu.
 - Click on **Device Customization**.
- PrecisionWeb Users:
 - Double click on the desktop shortcut or Internet Explorer Favorites for **i-STAT Customization**.


3.5: Update the CLEW and JAMS versions in the Customization Workspace.

- Under the “Default customization profile:” column, click on the “**i-STAT Analyzer CLEW**” button.

Note: Customization screens may vary depending upon i-STAT/DE version in use.

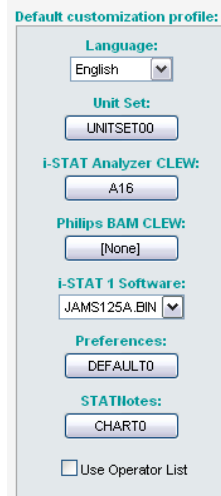

- Check the box next to the new version of **CLEW** and click **OK**.

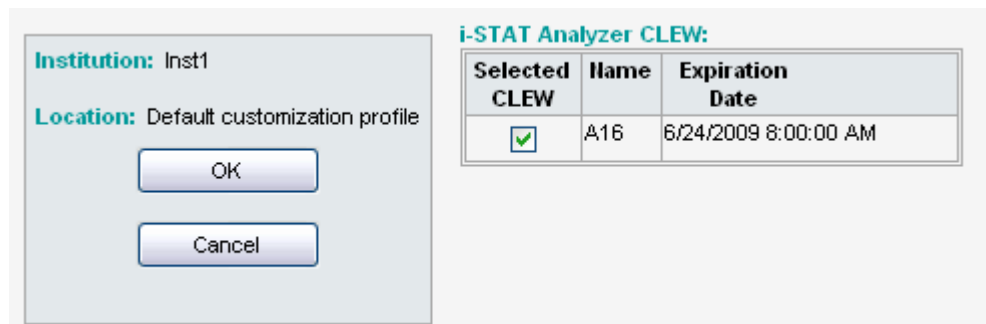

Answer **OK** to the question that appears.

- If “Uses Default” is not checked beside any Location-based customization profile, click the box under the “i-STAT Analyzer CLEW” column.

Location-based customization profiles:

| Location | Enabled | Uses Default | Update CLEW | i-STAT Analyzer CLEW | Philips BAM CLEW | Preferences | STATNotes | |
|----------|--------------------------|-------------------------------------|-------------------------------------|----------------------|------------------|-------------|-----------|--------------------------|
| ER | <input type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | A16 | [None] | DEFAULTO | CHARTO | <input type="checkbox"/> |
| ICU | <input type="checkbox"/> | <input type="checkbox"/> | <input checked="" type="checkbox"/> | A16 | [None] | DEFAULTO | CHARTO | <input type="checkbox"/> |
| Lab | <input type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | A16 | [None] | DEFAULTO | CHARTO | <input type="checkbox"/> |
| OR | <input type="checkbox"/> | <input type="checkbox"/> | <input checked="" type="checkbox"/> | A16 | [None] | DEFAULTO | CHARTO | <input type="checkbox"/> |

- Click the new version of **CLEW**, and then click **OK**.

Institution: Inst1

Location: Default customization profile

i-STAT Analyzer CLEW:

| Selected CLEW | Name | Expiration Date |
|-------------------------------------|------|----------------------|
| <input checked="" type="checkbox"/> | A16 | 6/24/2009 8:00:00 AM |

Answer OK to the question that appears.

- Under the “Default Customization profile:” column, click on the **i-STAT 1 Software** drop-down list. Select the JAMS version that matches the Product Update and click **OK**.

Default customization profile:

Language:
English

Unit Set:
UNITSET00

i-STAT Analyzer CLEW:
A16

Philips BAM CLEW:
[None]

i-STAT 1 Software:
JAMS125A.BIN

Preferences:
DEFAULTO

STATNotes:
CHARTO

Use Operator List

3.6: Enable Customization.

- If the **Enable Customization** box is not already checked, click the box next to this listing.
- Under the “Location-based customization profile:” section, make sure **Enabled** is checked for every location from which you wish to perform software updates on your i-STAT 1 Analyzers.

3.7: Update the software in the i-STAT 1 Analyzer.

- Go to the location where the i-STAT 1 Analyzer(s) you wish to update are located or contact someone at that location who can assist in updating the analyzer(s).

- Press the **On/Off** button on the analyzer.
- Press the **Menu** key to bring up the Administration Menu.
- Press **7 – Utility**. When prompted for a password, press **ENT**. If that did not work, a password is needed. Enter the password defined by your facility and press **ENT**. Note: Abbott Point of Care Inc. recommends changing the default password.
- From the Utility menu, press **3 – Receive Software**. A “Waiting to Send” message will appear on the analyzer display.
- Place the analyzer in the downloader or downloader / recharger. **Do NOT move the analyzer until step 3.8**. A **Communication in Progress** message will appear on the screen. After this disappears, the analyzer display will stay blank for approximately 5-10 seconds.
- The analyzer will then display 1’s and 0’s streaming across the screen signifying that it is receiving the software. Once the 1’s and 0’s disappear, the analyzer display will again go blank for approximately 5-10 seconds.
- A **Waiting to Send** message following by a **Communication in Progress** message will then appear on the analyzer display. After these messages disappear, the analyzer display will go blank, and the update process is complete.

3.8: Run the Electronic Simulator in the analyzer. When the simulator finishes, **PASS** should be displayed.

Note: If **PASS** is not displayed, re-run the Electronic Simulator. If the repeated Electronic Simulator attempt fails, contact APOC Technical Support. For additional information on running the Electronic Simulator, see Section 14 of the i-STAT 1 System Manual.

Congratulations. The process for updating the first i-STAT 1 Analyzer is complete.

Review the options below for additional instructions.

- If there are no additional analyzers to update, the process is complete.
- If there are additional analyzers to update, return to step **3.7**.

© 2022 Abbott. All rights reserved. All trademarks referenced are trademarks of either the Abbott group of companies or their respective owners.

This page intentionally left blank



i-STAT 1 SYSTEM TECHNICAL BULLETIN

Support Services

Abbott Point of Care and its distributors are committed to helping you resolve any problems with the i-STAT System: i-STAT 1 Handheld, cartridges, accessories and data management software. For technical assistance within the United States, please call Technical Services at 800-366-8020 toll free. Outside the U.S., please contact your local i-STAT distributor.

NORTH AMERICA

USA

Abbott Point of Care

400 College Road East
Princeton, NJ USA 08540
Tel: +800-366-8020, option 1
Tel: +800-284-0702 (waived customers)
Email: techsvc@apoc.abbott.com

CANADA

Abbott Point of Care
400 College Road East
Princeton, NJ USA 08540
Tel: +800-366-8020, option 1
Email: techsvc@apoc.abbott.com

AFRICA

ANGOLA

Promed International Health Ltd.
Largo AmilcarCabral n 2,A-B-C e 12 R/C
Ingombotas
Luanda, Angola
Tel.No.: +244 94 6146050

BOTSWANA

Leading Edge Markets (Pty) Ltd
P.O. Box 40551
Plot 128 Unit 5
Gaborone International
Finance Park
Gaborone, Botswana

DEMOCRATIC REPUBLIC OF CONGO Wagenia

259 Avenue Wagenia
Kinshasa 02476
Democratic Republic of Congo
Tel.No.: +243 9 97 000097

EGYPT

**International Company for Medical
Equipment S.A.E.**
24th, Gamiet El Dewal El Arabia
Giza, Egypt
Tel.No.: +202 33 47 06 78

KENYA

Phillips Healthcare Technologies Ltd
Phillips Business Park
Mombasa Road
Nairobi, Kenya
Tel.No.: 254 7336 12025

KENYA

Ziwala Limited
P.O. Box 17919
Nairobi, Kenya
00100
Tel. No: +254 722 725 529

LESOTHO AND MALAWI

Obsidian Health Ltd
 Cosmo Business Park
 Malibongwe Drive
 Randburg, South Africa 2188
 Tel.No.: +27 87 3535600

MAYOTTE AND REUNION ISLAND

Abbott France S.A
 40/48 rue d'Arcueil
 94593 Rungis Cedex, France
 Tel.No.: +33 1 4560 2500

NIGERIA

JB Consulting (MDP) Ltd
QDT Solution
 111 Concorde Road, Building 3,
 Level 1 Regus
 Maidenhead, Berkshire
 SL6 4BY UK
 Tel.No.: +44 7808 589 217

SUDAN

Penotee Multi Activities Co Ltd.
 Building No. 556
 Block 22, Badr Street
 Altaif, Khartoum Sudan
 Tel. No: +249 183 287794

ZAMBIA

Phillips Pharmaceuticals (Zambia)
 Plot #7236, Njolwe Road
 Light Industrial Area
 Lusaka, P.O. Box 36017
 Zambia

LIBYA

Al- Haremeen Pharmaceutical
Medical Supplies
 Extension of Al Saraj
 Street After Audi Showroom
 Tripoli, Libya
 Tel.No.: 201 00 3333 444

MOROCCO

Masterlab Sarl
 22 Rue IBN Tayeb Kadiri Les
 Orangers Rabat
 Morocco
 Tel. No.: +212 661208038

RWANDA

Phillips Pharmaceuticals (Rwanda)
Ltd
 GF 68, Kigali Modern City Market
 Commercial Street
 Kigali, Rwanda

TANZANIA

Phillips Pharmaceuticals (Tanzania)
Ltd
 Vingunguti Industrial Area
 P.O. Box 737
 Dae Es Salam, Tanzania
 Tel. No.: + 255 782637336

ZIMBABWE

Healthyard Laboratories Pvt Ltd
 Cooksey House
 2 Simon Mazorodze Rd., Southerton
 Harare AZ1095
 Zimbabwe
 Tel. No.: 263778401114

MAURITIUS

Health Focus Ltd.
 33 St. Clement Street
 Curepipe 74208
 Mauritius
 Tel.No.: +230 674 1000

NIGERIA

Phillips Pharmaceuticals
(Nigeria) Ltd
 122-132 Afprint Industrial Estate
 Apapa Oshodi Expressway
 Iyana-Isolo, Nigeria
 Tel.No.: +234 8056292422

SOUTH AFRICA

Obsidian Health Ltd
 Cosmo Business Park
 Malibongwe Drive
 Randburg, South Africa 2188
 Tel.No.: +27 87 3535600

UGANDA

Phillips Pharmaceuticals
(Uganda) Ltd
 N-S Maghavani Complex
 5th Street, Industrial Area
 Bugolobi
 Kampala, Uganda
 Tel.No.: 256 772 518 681

ASIA/PACIFIC (EXCLUDING CHINA)

AUSTRALIA

Abbott Australia
 299 Lane Cove Road
 Macquarie Park
 NSW 2113
 Australia
 Tel: 61 2 9857 1111

AUSTRALIA

CMM Technology
 Unit 12, 3 King Edward Road
 Osbourne Park
 Western Australia 6017
 Tel: +61 1300 79 70 30

BANGLADESH

Unimed Limited
 Rangs Nasim Square (7th Floor),
 House 46
 Sheikh Kamal Saroni, Road No. 16
 Dhanmondi
 Dhaka 1209 Bangladesh
 Tel: +8802 9128192-3

BRUNEI

Transmedic Pte Ltd.
 5, Jalan Kilang Barat
 9th Floor Petro Centre
 Singapore 159349
 Tel: +65 6 7371 945

COOK ISLANDS, FIJI, FR. POLYNESIA, SAMOA, TONGA, AND TOKELAU

Abbott New Zealand
 Building D, 4 Pacific Rise
 Mt. Wellington Auckland 1060
 New Zealand
 Tel: 64 9 573 6030

HONG KONG

BioAsia Diagnostics Company Ltd
 Unit 1-3, 20F
 CRE Centre, 889 Cheung
 Sha Wan Road, Cheung Sha Wan
 Kowloon Hong Kong
 Tel: +852 27870906

INDIA

Sandor Medicaids Private Ltd.
8-2-326/5, Road No. 3
Banjara Hills
Hyderabad, India 500 034
Tel: +91 40 233 570 48

INDONESIA

Transmedic Pte Ltd.
5, Jalan Kilang Barat
9th Floor Petro Centre
Singapore 159349
Tel: +65 6 7371 945

JAPAN

Japan Medical Link Co., LTD
8-1 marunouchi, matsumoto-shi
nagano-ken, JP 390-0873
Tel: +81 263 38 0411

JAPAN

Konishi Medical Instruments Co., Ltd
2-1-5awajimachi, chuo-ku
osaka, JP 541-0047
Tel: +81 6 6941 1363

JAPAN

Masuda Medical Instruments Co., LTD
50 takedawarayacho, hushimi-ku kyoto,
JP 612-8443
Tel: +81 75 623 7111

JAPAN

Nakakita Yakuin Corporation
3-5-15 marunouchi
naka-ku
aichi-ken, JP 460-8515
Tel: +81 52 971 3681

JAPAN

Takeyama Co., LTD
16-1-5 kita6jonishi
chuo-ku
hokkaido, JP 060-0006
Tel: +81 11-0611-0100

JAPAN

Abbott Japan
3-5-27, Mita, Minato-ku
Chiba 270-2214
Japan
Tel: 65 6914-8351

INDIA

Crescent Organics Pvt. Ltd.
WINDSOR, 2nd Floor
Cst Road, Kalina, SantaCruz (E)
Mumbai 400098

JAPAN

Hokuyaku Inc.
16-1-5 kita6jonishi, chuo-ku
Hokkaido JP 060-0006
Tel: +8 11611-0989

JAPAN

Kawanishi Corporation
1-3-9 ima, kita-ku
okayama-ken, JP 700-8577
Tel: +81 86 241 9242

JAPAN

Kyowa Medical Corporation
156-2ikedada, suruga-ku
shizuoka-ken, JP 422-8005
Tel: +81 45 655 6600

JAPAN

Mitas, Inc.
4-901 tonyacho
hukui-shi
hukui-ken, JP 918-8556
Tel: +81 776 28 2888

JAPAN

Sanshodo Co., LTD
68kamitobaomonocho
minami-ku
kyoto, JP 601-8533
Tel: +81 75 681 5131

JAPAN

Tomiki Medical Instruments, Co., Ltd
2-46 tonyacho
kanazawa-oshi
Ishikawa-kens, JP 920-8539
Tel: +81 76 237 5555

**KIRIBATI, MICRONESIA,
NEW CALEDONIA, AND
PAPUA NEW GUINEA**

Abbott Australia
299 Lane Cove Road
Macquarie Park
NSW 2113 Australia
Tel: 61 2 9857 1111

INDIA

Genpharmasec Limited
104/105, 1st Floor
Gundecha Indust. Comp. Premise
Akurli Road, Kandivali East
Mumbai II 400 India

JAPAN

Iwabuchi Pharmaceutical Co., Ltd
1-5 Takanodai
Yotsukaido-Shi (Chiba) JP
284-0033
Tel: +81 43 310 6630

JAPAN

Kishiya Inc
1-14-21 matsushima
higashi-ku
hukuoka-ken, JP 812-0062
Tel: +81 92611 8000

JAPAN

Maruki Medical Systems Inc.
3-20-7 nishinakata, taihaku-ku
miyagi-ken, JP 981-1105
Tel: +81 22 242 6001

JAPAN

Musashi Engineering Co., LTD
sekitori60, uedacho
Aichi, JP 441-8134
Tel: +81 532 25 1231

JAPAN

Shibataintech
2-11-3 oroshimachi
wakabayashi-ku
miyagi-ken, JP 984-0015
Tel: +81 22 236 2311

JAPAN

**Yamashita Medical Instruments
Co., LTD**
hakataza and nishigin building
10
2-1 shimokawabatacho
fukuokaken, JP 812-0027
Tel: +81 956 25 2112

KYRGYZSTAN

Meditek Ltd.
Ryskulova Str. 12
Bishkek, Kyrgyzstan 720001
Tel.: 00996 312 93 50 05

MALAYSIA

Transmedic Pte Ltd.
5, Jalan Kilang Barat
9th Floor Petro Centre
Singapore 159349
Tel: +65 6 7371 945

NEPAL

Synergy Medisales Pvt. Ltd
#205, Bishal Nagar, Bishal Basti kha
Kathmandu, B
977 Nepal
Tel.: +977-01-4425956

PAKISTAN

Hospicare Systems
Mezzanine Floor. Rabbiya Garden Block 3
Karachi
74800, Pakistan
Tel. No: +92 2134921913

SINGAPORE

Transmedic Pte Ltd.
5, Jalan Kilang Barat
9th Floor Petro Centre
Singapore 159349
Tel: +65 6 7371 945

SRI LANKA

Hemas Surgicals & Diagnostics Pvt. Ltd
No. 75 Braybrook Place
Colombo, Sri Lanka
Tel: 94114766680

THAILAND

Connect Diagnostics Co., Ltd
1/57-58 Venice Di Iris
Soi Watcharaphol 2/7
Bangkok 10220 Thailand
Tel: 6623470109

MALDIVES

Medtech Maldives Pvt Ltd
Ma. Fas Eri 3rd Floor
Ameenee Magu Male 20205
Male 20205 Maldives
Tel: +960 3308663

NEW CALEDONIA

Medi-Services S.A.R.L
8, rue Reverce
Noumea
New Caledonia
Tel: 687272000

PHILIPPINES

Euromed Laboratories Phil., Inc.
Sabutan Paliparan Road
Brgy. Sabutan
Silang. Cavite 4118 Philippines
Tel.: 63 (02) 85240091 98

**SOLOMON ISLANDS,
TUVALU, VANUATU, AND YAP**

Abbott Australia
299 Lane Cove Road
Macquarie Park
NSW 2113
Australia
Tel: 61 2 9857 1111

TAIWAN

Ditech Enterprise Co., Ltd.
9F, No. 168, Sec. 1
Zhongshan Rd
New Taipei City 234 – TA - 234
Taiwan
Tel: +886 2 8923 2236

UZBEKISTAN

Albatros World
Buyuk Ipak Yoli street, 49
Mirzo Ulugbek district
Tashkent, Uzbekistan 10005 4
Tel: +998 97 744 78 63

MYANMAR (BURMA)

**Advanced Diagnostic Products
Trading Co., Ltd**
No. 34, Thate Pan Street
Ahlone Township
Yangon, Myanmar
Tel.: +951 229587

NEW ZEALAND

Abbott New Zealand
Building D, 4 Pacific Rise
Mt. Wellington
Auckland 1060
New Zealand
Tel: 64 9 573 6030

PHILIPPINES

**Medical Trends & Technologies,
Inc. (MEDT2EK)**
#200J Elizalde Street, Phase V
BF Homes
Paranaque City 1700
Philippines
Tel: 6388069767

SOUTH KOREA

POC Korea Ltd.
Ace Gasan Tower, 121 Digital-ro
Suite 1501
Gumchon-gu, Seoul
08505 South Korea
Tel: +82 2 2058 0971

THAILAND

Transmedic Pte Ltd.
5, Jalan Kilang Barat
9th Floor Petro Centre
Singapore 159349
Tel: +65 6 7371 945

VIETNAM

Transmedic Pte Ltd.
5, Jalan Kilang Barat
9th Floor Petro Centre
Singapore 159349
Tel: +65 6 7371 945

CHINA

BIOTEC Biotechnology Limited Co.
Room 1013, Wu Yang Xin Chen Plaza
No. 111-115, Si You Xing Road
Yue Xiu District, Guangzhou
People's Republic of China
Tel. No: 86-020-28865008

**Beijing Huiwen Yuanmei
Technology & Trading Co., Ltd**
Room C100, Area C, No. 14 Building
No. 2 Kechuang Dongwu Street
Beijing, China 100111
Tel. No: 010-57621837

**Guangdong Biochem
Healthcare Investment Co., Ltd**
Room 903-904, Building A, No. 263
Shuangcheng International Plaza
Guangzhou, China
Tel. No: 86-020-38372808

**Chongqing Chenjie Medical
Equipment Co. Ltd**
Room 4-6, Building 4
No. 101 Cuibai Rd
Dadukou District
Chongqing, 400080 China
Tel.No: 86-023-68153822

**Sinopharm Group Hebei
Medical Equipment Co., Ltd.**
No. 309 Zhonghua North Street
Xinhua District
Shijiazhuang 050050 China
Tel. No: 0311-67302513

**Sinopharm Medical Supply
Chain Service (AnHui) Co., Ltd**
No. 15, Red Maple Road
Hi-Tech Zone
HeFei
230088 China
Tel. No: 86-551-65310478

**Sinopharm Group Xinjiang Medical
Equipment Co., Ltd**
16th Floor, 1st Building
199 North Liyushan Road
Urumai 830054 China
Tel. No: 86-991-3197595

**Xuzhou Ruida Medical Equipment
Ltd**
Room 306, No. 25, Jinqiao Road
Jinshan Bridge Economic
Development Zone
Xuzhou 221000 China
Tel. No: 86-18761439568

**Shenzhen Zichang Technology
Co., Ltd.**
Room 603, 6 Building, Tongfuyu
industrial city
Tanglang, Xili Street, Nanshan
District
Shenzhen CH
Tel. No: 18124785081

**Beijing Bao Ze Kang
Pharmaceutical Co., Ltd.**
B512/513 5F, 6-1 building
28 red lotus road
Xicheng District
Beijing CH
Tel. No: 010-62358509

**Chengdu Life Medical Equipment Co.
Ltd.**
Room 2-5, 13th Floor, Building No.2
799 Jinfu Rd
Jinniu District, Chengdu city
Sichuan Province
People's Republic of China
Tel. No: 028-65037108

**Realcan Phamaceutical Group
Shanghai Co., Ltd.**
Room 1501, 388 Fenglin Road
Xuhui District
Shanghai
Shanghai CH 200032
People's Republic of China
Tel. No: 86-18153573207

**Wuhan Zhongji Pathology Diagnosis
Center Co., Ltd.**
Room 1, 4th Floor, Building B12
No. 818 High-tech Street
EastLake Development Zone
Wuhan, China
Tel. No: 15071168946

**Shanghai Jointown Medical Devices
Supply Chain CO.,LTD**
Room 206, 2/F, Building 7
No. 285 fengguan Road
Jinshan District
Shanghai, CH
Tel. No: 15207181724

**Guangdong heyidukang Medical
Technology Co., Ltd.**
Room 217, building 7, Jiangmen
Wanda Plaza, Pengjiang district,
Jiangmen
China 529000
Tel. No: 15807502929

**Sinochem Beijing Medical Technology
Co., Ltd**
17 /F, Jincheng Building
Longpan Middle Road
Qinhuai District, Nanjing, China

EUROPE

AUSTRIA

Abbott Austria
Perfektastr., 84A
Wien, Austria
A 1230
Tel. No.: +43 1 89 122 0

AZERBAIJAN

Albatros Healthcare
Y. Haseynov 6
Narimanov
Baku, Azerbaijan AZ1021
Tel.No.: +994 125648635

BELGIUM & LUXEMBOURG

**Abbott sa/nv
(Abbott Belgium)**
Avenue Einstein 14
1300 Wavre, Belgium
Tel.No.: +32 10 47 53 11

CROATIA AND MALTA

Mark Medical d.o.o.
Budmanijeva 5
Zagreb, 10000, Croatia
Tel. No.: +385 1 6065 444

CZECH REPUBLIC

Coro Medical s.r.o.
Biskupsky dvur 2095/8
Prague, Czech Republic 110 00
Tel.No.: 42603229000

DENMARK, GREENLAND, ICELAND AND FAROE ISLANDS

**Abbott Laboratories A/S
(Abbott Denmark)**
Emdrupvej 28C
DK-2100
Copenhagen, Denmark
Tel. No.: +45 39 77 00 00

FINLAND

Abbott Oy (Abbott Finland)
Linnoitustie 4
Espoo, Finland 02600
Tel. No.: +358 9 751 8418

FRANCE

Abbott France S.A
40/48 rue d'Arcueil
94593 Rungis Cedex, France
Tel. No.: +33 1 4560 2500

GEORGIA

**Advanced Medical Technologies &
Services AMTS**
Saburtalo Str. 7, Flat 13
Tbilisi 0159
Georgian Republic
Tel.No.: +995 32 532 032/036

GERMANY

Abbott GmbH (Abbott Germany)
Max-Planck-Ring 2-3
65205 Wiesbaden – Delkenheim
Germany
Tel.No.: +49-6122-58-0

GERMANY

WS Laborbetriebsgesellschaft GmbH
Langendembach 72A
Langenoria GE 07381
Germany
Tel.No.: +49 3647 – 414378

GREECE, CYPRUS, ALBANIA AND KOSOVO

Pegasus Health and Sciences S.A.
Thessalonikis 18
Chalandri
Athens, GR 15234
Tel.No.: +30 210 6139496

IRELAND

Abbott Laboratories, Ireland Limited
Block B, Liffey Valley Office Campus
Dublin 22
D22 XOY3
Tel. No.: +353 1 469 1560

ITALY, VATICAN CITY, MALTA AND SAN MARINO

Abbott SRL (Italy)
Via Ribotta 9
Rome, Italy 00144
Tel. No.: +39 06 52 9911

ITALY

Burke & Burke
Via A Einstein, 32
Assago (MI)
20090 Italy
Tel. No.: +39 02 4571 3663

ITALY

Tobar Service SRL
Via Veronica Gambarà 23
Rome, Italy
Tel. No.: +39 06 89161897

NETHERLANDS

**Abbott B.V. (Abbott The
Netherlands)**
Wegalaan 9
2132 Hofddorp
Netherlands
Tel.No.: +31 88 82 22 500

NORWAY

**Abbott Norge as (Abbott
Norway)**
Martin Linges vei 25
Postboks 1 Fornebu N-1330
Norway
Tel. No.: +47 81 55 99 20

POLAND

Diag-med
ul.Stanisława 50
Prusków
05-800 Poland
Tel.No.: +48 22 8389723

PORTUGAL AND SPAIN

Magnamed, Lda
Rua Prof Francisco Gentil n° 22F
2620-097 Póvoa de Santo Adrião
Portugal
Tel.No.: +351 21 938 32 40

ROMANIA

**MD&D (SC Medical Devices &
Diagnostics SRL)**
111C, Drumul Vales Furcii St.
Sector 6
Bucharest 061985
Romania
Tel.No.: 40311054280

ROMANIA

Ducos Trading Ltd
Johann Strauss, No. 2A
Etaj 3, Camera 5A
Bucuresti 020312
Romania
Tel. No: +40 314 381 094

RUSSIA

Eco-med-s M LLC
Samokatnaya str., 2A
Bld. 1 p.105, r. 3, of. 2
Moscow, Russia 111033
Tel. No: +7 (495) 748-43-50

**SPAIN, PORTUGAL, GIBRALTAR,
ANDORRA & CANARY ISLANDS**

**Abbott laboratories SRL
(Abbott Spain)**
Parque Empresarial Via Norte
c/ Quintanavides
17 Madrid 28050 Spain
Tel. No.: +34 91 337 3400

**SWITZERLAND AND
LIECHTENSTEIN**

Axonlab AG - Switzerland
Taefernstrasse 15
CH-5405 Baden-Daetwil
Switzerland
Tel.No.: +41 56 484 8080

UNITED KINGDOM

**Abbott Laboratories Limited
(Abbott UK)**
Abbott House
Vanwall Business Park
Maidenhead
SL6 4XF, United Kingdom
Tel.: 44 16 28 784 041

RUSSIA

**Abbott Laboratories Ltd (Abbott
Russia)**
Leningradskoe shosse 16A
Bld 1, Moscow
Russia 125171

RUSSIA

UDS (Updated Diagnostic Systems)
3 Bogorodsky Vall
Moscow, Russia 107076

SPAIN

Technédica
C/ Luis Doreste Silva, 62 - 1°
Las Palmas
35004 Spain
Tel. No.: +34 928 242381

**SWITZERLAND AND
LIECHTENSTEIN**

Abbott AG (Abbott Switzerland)
Neuhofstrasse 23
CH-6341 Baar
Switzerland
Tel. No: +41 41 768 44 44

RUSSIA

Servis Instrument LLC
Bld 1, Room 18
Ryabinovaya Str. 38
Moscow 121471
Russia
Tel. No: +7 495 628 7845

**SLOVENIA, BOSNIA, AND
HERZEGOVINA**

EUROMED, d.o.o.
Podpeska cesta 14
Brezovica pri Ljubljani
Slovenia, SI-1351
Tel. No: +386 59 338 202
+386 40 451 115

SWEDEN

**Abbott Scandinavia AB (Abbott
Sweden)**
Hemvärnsgatan 9
Solna SE17129
Sweden
Tel. No.: +46 (0)8 546 567 18
E.post: Service@abbott.se

TURKEY

**Gazi Kimya Tip Teknolojileri
San. Ve Tic. Anonim Sirketi**
Serifali Mah. Beyan Sokak
Akalin Plaza, No. 3
Umraniye, Istanbul Turkey
Tel. No.: 0090 532 462 55 69

LATIN AMERICA

ANGUILLA

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel. No: 407 475 1168

ARUBA

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

BELIZE

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

BONAIRE

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

BRAZIL

Medcare Produtos
E Equipamentos Medico Hospitalares LTDA
Al Tres de Outubro, 795
Sarandi-Porto Alegre
RS, Porto Alegre
91130-470, Brazil
Tel. No: +51-984466491

CHILE

Alatheia Medical SpA
Avenida del Valle Norte 945
Oficina 5610
Huechuraba, Santiago
Chile
Tel.: +562 343 5122

COLOMBIA

Comprolab
CALLE 106 #54-63
Bogota
Colombia

DOMINICA

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

ANTIGUA AND BARBUDA

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

THE BAHAMAS

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

BERMUDA

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

BRAZIL

Tecno4 Produtos Hospitalares Ltda.
Rua Backer
89 Cambuci
Sao Paulo
Brazil 01541-000
Tel.: +1 55 11 95557-6345

BRITISH VIRGIN ISLANDS

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

COLOMBIA

Abbott Laboratories de Colombia S.A.S
Calle 100 #9A-45
Piso 14 Bogota
Colombia
Tel.: (57 1) 628 5600

COSTA RICA

Tri Dm
Calles 16-18 Ave 8
No. 1628
San Jose, CR 10103
Tel.: +506 2257 7676

DOMINICAN REPUBLIC

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

ARGENTINA

Drogueria Artigas S.A.
Av. Jose Luis Chorroarin, 1079
Buenos Aires Capital Federal
Argentina C1427CXH
Tel.: +64 3 338 0999

BARBADOS

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

BOLIVIA

Sanimed S.R.L.
Calle Cordero No 130 Zona
San Jorge, LaPaz, Bolivia
Tel.: 591-2 2431294

BRAZIL

PMH Produtos Medicos
Hospitalares LTDA
SIA Trecho
17 RUA 8 LOTE 170
Brasilia, DF Brazil 71200-222
Tel.: +55 61 81312112

CAYMAN ISLANDS

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

COLOMBIA

Arrow Medical SAS
Calle 23# 43 A-100
Medellin, Antioquia
Colombia
Tel.: +574 356 1111

CURACAO

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

ECUADOR

M. R. U.
Av. Eloy Alfaro N29-235
Entre Italia y Alemania P
Quito, Ecuador
Tel.: +593-2 3801 533

EL SALVADOR**Servicios Quirúrgicos**

Km. 11 Carretera al Puerto de la Libertad
 Antigua Cuascatlan
 La Libertad, El Salvador
 Tel.: +503 2241-6200

EL SALVADOR**Interserv SA de CV - El Salvador**

67 Avenida Sur
 Final Pasaje "A", No. 4-B
 Colonia Roma
 El Salvador Cp 1101
 Tel: +503 2223-5510

GRENADA**American Hospital Supply (AHS)**

1060 Maitland Center Commons
 Maitland, FL USA 32751
 Tel.: 407 475 1168

GADELOUPE, GUYANA, AND MARTINIQUE**Abbott France S.A**

40/48 rue d'Arcueil
 94593 Rungis Cedex, France
 Tel.No.: +33 1 4560 2500

GUATEMALA**Promed S.A.**

Avenida Ferrocarril 19-97
 Empresarial El Cortijo Interior 202
 Zona 12, Guatemala 10001
 Tel.: +502 2223 86 00

GUYANA**American Hospital Supply (AHS)**

1060 Maitland Center Commons
 Maitland, FL USA 32751
 Tel.: 407 475 1168

HAITI**American Hospital Supply (AHS)**

1060 Maitland Center Commons
 Maitland, FL USA 32751
 Tel.: 407 475 1168

HONDURAS**ProdyLab S de R.L.**

Col. Santa Ana, 11 y 12
 Avenida 12 Calle Bloque No. 117,
 N.O Sand Pedro Sula, Honduras
 Tel.: +504 2550-3091

JAMAICA**American Hospital Supply (AHS)**

1060 Maitland Center Commons
 Maitland, FL USA 32751
 Tel.: 407 475 1168

MEXICO**Abbott Laboratories de Mexico SA de CV**

Calzada de Tlalpan, 3092
 Colonia Ex Hacienda Coapa
 D.F. 04980 Mexico
 Tel.: (5255) 58097500

MEXICO**Abalat S.A. de C.V.**

Abasolo 78
 Pedregal de Sta Ursula
 Coyoacán, Mexico City MX 04650
 Tel.: 80001500

MEXICO**Casa Plarre S.A. de C.V.**

Av. Cuauhtemoc 220, interior
 201 Colonia Doctores
 Alcaldía Cuauhtemoc
 Ciudad de Mexico

MEXICO**Centrum Promotora Internacional S.A. de C.V.**

Medellin 324
 Roma Sur, Cuauhtémoc
 Mexico City, MX 06760
 Tel.: +52 5552652500

MEXICO**Comercializadora de Reactivos**

Ave. Eugenio Garza Sada N°.6336
 Col. Ciudad Satélite
 Monterrey, MX 64960
 Tel.: 800 7107 566

MEXICO**Diagmex S.A. de C.V.**

Calle San Martin de Porres 3777
 y 3769
 Jardín de San Ignacio
 Zapopan, Jalisco MX 45040
 Tel.: 01 8008770007

MEXICO**Diagnóstica Internacional**

Circunvalación Sur 156
 Las Fuentes
 Zapopan, Jalisco MX 45070
 Tel.: +52 3337771940

MEXICO**GRUPO ICTZE SA DE CV**

Callejon SN 172 Parcela 29
 Col. Playa de vacas
 Medellin, Veracruz MX 94274
 Tel.: 2299357143

MEXICO**Grupo Biomédico Empresarial**

Calle Tlacoquemecatli 21
 Tlacoquemecatli del Valle
 Benito Juarez
 Mexico City MX 03200

MEXICO**Impromed S.A. de C.V.**

Perpetua 30, Int. 101
 Col san José Insurgentes
 Alcaldía Benito Juarez
 Mexico City MX 03900
 Tel.No: +52 5553228870

MEXICO**Instrumentos y Equipos Falcón S.A de C.V**

Calle José María de Teresa 188
 Campestre – Álvaro Obregón
 Mexico City MX 01040
 Tel.No: 52 55 17190600

MEXICO**LAVKEM SAPI de CV**

Paseo de los Tamarindos, 90
 Torre 1, Piso 32
 Ciudad de Mexico

MEXICO

Ruvel S.A. de C.V.
Av Lago de Guadalupe KM 1.5
San José el Jaral
Cd López Mateos, MX 54927
Tel.No: +52 5553701047

PANAMA AND GUATEMALA

Promed S.A.
Pque. Industrial Costa del Este,
C/ 2da Edi.Premed
Panama
Tel.: +507 303-3100

PUERTO RICO

Cardinal Health Puerto Rico
Centro Internacional de Distribution
Carretera 869 KM 4.2
Guaynabo 00695 Puerto Rico
Tel.: 787 625 4100

ST. MAARTEN

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

TORTOLA

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

URUGUAY

Bioerix S.A.
Simon Bolivar 1472
Montevideo Uruguay
Tel.: 59827081624

MEXICO

**SERVICIOS INTEGRALES E
INGENIERIA HENLIS**
Bosque de las Palmas
Naucalpan de Juárez, MX 52787
Tel.No: 52722 216 1905

PARAGUAY

Index S.A.C.I.
Boqueron No. 676
Postal Code 1414
Asuncion, Paraguay
Tel.: 595 21 214 213

ST. KITTS AND NEVIS

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

**ST. VINCENT AND THE
GRENADINES**

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

TRINIDAD AND TOBAGO

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

MONTSERRAT

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

PERU

Biomedical Systems, S.A.
Av. Rafael Escardo
No. 823, San Miguel
Lima 32 , Peru
Tel: +511 241 3843

ST. LUCIA

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

SURINAME

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

**TURKS & CAICOS
ISLANDS**

American Hospital Supply (AHS)
1060 Maitland Center Commons
Maitland, FL USA 32751
Tel.: 407 475 1168

MIDDLE EAST

IRAN

Paya Zist Arayeh
No. 13, 10th Alley
Zarafshan St., Farahzadi St.
Shahrak-e-Gharb, Tehran
Iran
Tel: 989122836214

KINGDOM OF BAHRAIN

Gulf Corporation for Technology
Bldg 2038, Road: 4156, Block: 341
Manama/Al Juffair
Kingdom of Bahrain
Tel: 97317240597

OMAN

Waleed Pharmacy & Stores LLC
Street No. 40, Building No. 176
North Al Ghoubra, Ministry Road
Ghoubra
Muscat 100 Oman
Tel: 96899338480

SAUDI ARABIA

**Medical Supplies & Services Co. Ltd
(MEDISERV)**
Mediserv Building
Alahassa Street
Riyadh, Saudi Arabia
Tel.No.: +966 1 478 0555

ISRAEL

Rhenium Ltd.
20 Hasatat st PO Box 180
Modi'in 71711 Israel
Tel: +972-8-9558888

KUWAIT

Bader Sultan & Bros. W.L.L.
Block B 67, Shuwaikh, Future Zone
13009 Safat
Kuwait
Tel: +965 2461 0280

PALESTINE

**Medical Supplies and
Services (MSS)**
Betunia, Unipal building
Ramallah Palestine
Tel.No.: 00972 2 2959372

U.A.E.

Gulf & World Traders
P.O. Box 5527, Al Garhoud Area
Dubai
Tel.No.: +971 4 2821717

JORDAN

Smart Care Medical Co.
Wasfi El tall street
al Fannar complex
Amman, Jordan
Tel: 1 55 11 95557-6345

LEBANON

Union Pharmaceutique d Orient
Museum Street Badaro
Beirut 11-0709
Lebanon
Tel: +96 113 87815

QATAR

**Khalid Scientific
Company W.L.L.**
Airport St., Bldg. No. 402
Area No. 47
Doha
Qatar 4349
Tel: 974 432 5198

YEMEN

Griffin Ltd.
Near Linton R/A 13
Opposite Mineral Circles
Bearing Free St.
Jebel Ali, 732 UAE
Tel. No.: +967 1 441 907

