AFINION™ 2 USB Flash Drive

REF 1117519

- en
- da
- de
- el
- es
- fi
- fr
- it
- nl
- no
- pt
- ru
- sv

1117532 Rev. A 2026/04

Abbott Diagnostics Technologies AS
Kjelsåsveien 161
P.O. Box 6863 Rodelekk
NO-0504 Oslo, Norway
www.globalpointofcare.abbott

ISO 13485 certified company

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en AFINION™ 2 USB Flash Drive

Intended use

Afinion® USB Flash Drive should be used to upgrade the software on Afinion analyzers.

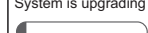
Kit contents

- 1 Afinion USB Flash Drive
- 1 Package Insert

Procedure

This procedure is valid for:

- **AFINION™ 2** Analyzer

<div>1</div> <div><div>Software version: X.XX New version: Y.YY Continue to upgrade? </div></div>	<ul style="list-style-type: none">• Switch off the Afinion analyzer• Insert the Afinion USB Flash Drive into one of the two available USB ports at the back of the analyzer.• Switch on the analyzer and follow the instructions on the touch screen.• To start the software upgrade, touch: • To abort the software upgrade, touch: 
<div>2</div> <div><div>X.XX System is upgrading </div></div>	<p>System is upgrading</p> <p>The software upgrade takes 4-5 minutes.</p>
<div>3</div> <div><div>Software upgrade completed. Remove flash drive to restart.</div></div>	<ul style="list-style-type: none">• When the SW upgrade is finished, remove the Afinion USB Flash Drive from the USB port at the back of the analyzer.• The analyzer will automatically restart and the new software version will be visible in the upper left corner of the display. The self-test will be performed. Do not turn the analyzer off before the self-test is finished.• The analyzer is now upgraded with the newest software version.

Troubleshooting

If the SW upgrade fails or the analyzer is giving **Self-test error - Analyzer in non-operative mode**, follow the instructions below:

- Switch off the analyzer by pulling out the plug from the power socket.
- Restart the analyzer and follow the SW upgrade procedure.
- If the SW upgrade fails again, contact your local Afinion supplier.

Consult the Afinion analyzer User Manual for detailed instructions on how to operate the analyzer.

Technical Support

Africa & RCIS (Russia and Commonwealth of Independent States	arcis.techsupport@abbott.com
APAC (Asia-Pacific)	AP.TechSupport@abbott.com
Europe & Middle East	EME.TechSupport@abbott.com
LATAM (Latin-America)	LA.TechSupport@abbott.com
United States	Afinion.support@abbott.com

For US customers only

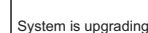
The manufacturer provides a toll free line for technical support.
Call: 1-866-216-9505. The toll free number is available for use only in the United States of America.

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© 2026 Abbott. All rights reserved. All trademarks referenced are trademarks of either the Abbott group of companies or their respective owners.